


Instructions for use

INTRAsurg 300 / INTRAsurg 300 plus



Always be on the safe side.

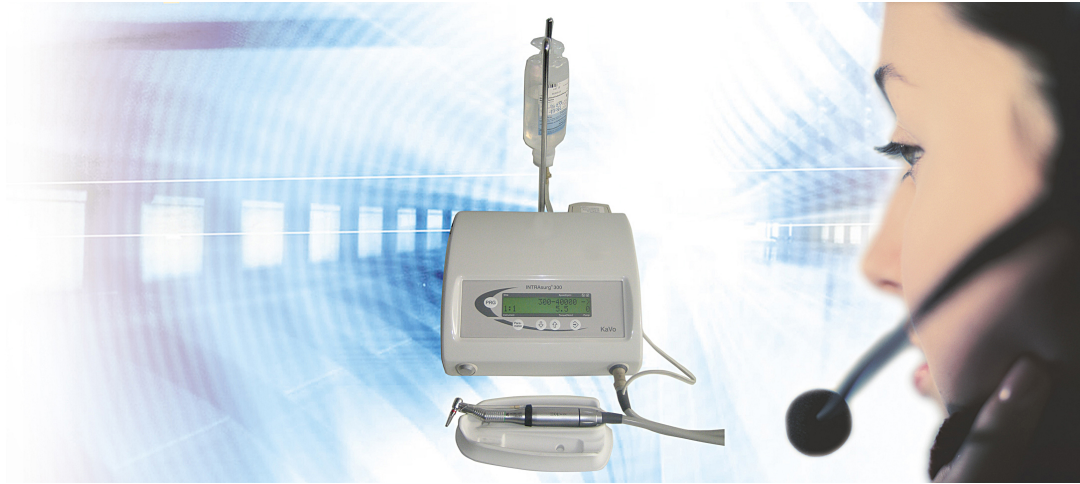

KaVo. Dental Excellence.

Distributed by:
KaVo Dental GmbH
Bismarckring 39
D-88400 Biberach
Tel. +49 7351 56-0
Fax +49 7351 56-1488

Manufacturer:
Kaltenbach & Voigt GmbH
Bismarckring 39
D-88400 Biberach
www.kavo.com



Contents

1	User instructions	3
1.1	User guide	3
1.1.1	Symbols	3
1.2	Target group	3
1.3	Service	3
1.4	Warranty terms and conditions	3
1.5	Transportation and storage	4
1.5.1	Currently valid packaging regulations	4
1.5.2	Damage in transit	4
1.5.3	Information on the packaging: Storage and transportation	5
2	Safety	7
2.1	Description of safety instructions	7
2.1.1	Warning symbol	7
2.1.2	Structure	7
2.1.3	Description of danger levels	7
2.2	Purpose – Proper use	7
2.2.1	Information on electromagnetic compatibility	9
2.2.2	Risks from electromagnetic fields (implanted life-sustaining system)	10
2.2.3	Disposal of electronic and electrical devices	10
2.3	Safety instructions	11
3	Product description	13
3.1	INTRAsurg 300/300 plus	13
3.2	Controls	14
3.3	Multifunctional foot control	15
3.4	Rating plate for 300/300 plus	16
3.5	Technical data for the 300/300 plus	16
4	First use	18
4.1	Electrical connection	18
4.2	Unpacking	19
4.3	Adjust the mains input voltage	19
4.4	Mount the bottle holder and use the correct coolant container	21
4.5	Connect the surgery motor and coolant reservoir	22
5	Operation	25
5.1	Operation in general	25
5.1.1	Start up	25
5.2	Surgical motors	25
5.2.1	INTRA S 550 surgical motor	25
5.2.2	INTRA LUX SL 550 surgical motor	26
5.2.3	Remove the handpiece or angle piece	27
5.2.4	Start the motor	27
5.3	Free use	28
5.3.1	Parameters menu	28
5.3.2	Set the torque limitation	30
5.3.3	Preset the speed range	30
5.3.4	Set the coolant flow	31
5.3.5	Selecting the handpiece and contra-angle handpiece	32
5.3.6	Switching the direction of motor rotation	33







Contents

5.3.7	Rinsing function	33
5.3.8	Surgical motor in free use	34
5.3.9	Request maximum torque	34
5.4	Program	35
5.4.1	Start the program	35
5.4.2	Switch forward step	36
5.4.3	Go backward a step	36
5.4.4	Set parameters	36
5.4.5	Limit the number of steps	37
5.4.6	Setting up delivery status	38
5.4.7	Surgical motor in the program	38
5.4.8	Request maximum torque	39
5.5	Functions exclusively with the 300 plus program version	40
5.5.1	Recognition function	40
5.5.2	Light functions	42
5.5.3	Request maximum torque (free use)	43
5.5.4	Request maximum torque (program)	43
5.5.5	Recommended programming when setting a number of implants one after the other	43
5.5.6	INTRAsurg calibration	44
6	Preparation methods DIN EN ISO 17664	46
6.1	Cleaning	47
6.1.1	Machine cleaning	47
6.1.2	Manual cleaning	47
6.1.3	Rinse the hoses	47
6.2	Disinfection	49
6.2.1	Automated disinfection	49
6.2.2	Manual disinfection	49
6.2.3	Drying	50
6.2.4	Service, inspection and testing after preparation	50
6.3	Sterilisation in a steriliser in compliance with DIN EN 13060	50
6.3.1	Packaging	51
6.3.2	Sterilisation	52
6.3.3	Storage	52
7	Troubleshooting	53
7.1	Exchange the pump hose	56
7.2	Exchange coolant hose	57
7.2.1	Exchange the coolant hose from the handpiece or angle piece.	58
7.2.2	Exchange the coolant hose in the motor hose	58
7.2.3	Disposable coolant hose	59
7.3	Device settings table	61
8	Accessories	62
9	Details on electromagnetic compatibility	63
9.1	Guidelines and manufacturer's declaration - electromagnetic transmission	63
9.2	Guidelines and manufacturer's declaration - electromagnetic resistance to jamming	63
9.3	Guidelines and manufacturer's declaration - electromagnetic resistance to jamming	64
9.4	Recommended safe distance between portable and mobile HF telecommunications equipment and the INTRAsurg 300	65

1 User instructions

1.1 User guide

1.1.1 Symbols

	Refer to chapter Safety/Warning
	Important information for users and technicians
	Thermodisinfected
	Sterilisable up to 135°C
	CE mark (Communauté Européenne). A product with this mark meets the requirements of the relevant EC directives, i.e. the applicable standards in Europe.
	ESD warning plate

1.2 Target group

This document is for dentists and office personnel.

1.3 Service



Service hotline:

+49 (0) 7351 56-2500

Service.Einrichtungen@kavo.com

Please indicate the product serial number in all requests.

Additional information can be obtained at: www.kavo.com Additional information can be obtained at: www.kavo.com

1.4 Warranty terms and conditions

KaVo provides the end customer with a warranty that the product cited in the hand-over certificate will function properly and guarantees zero defects in the material or processing for a period of 12 months from data of purchase, subject to the following conditions:

Upon justified complaints of flaws or a short delivery, KaVo will make good its warranty by replacing the product free of cost or repairing it according to the customer's wishes. Other claims of any nature whatsoever, in particular with respect to compensation, are excluded. In the event of default and gross negligence or intent, this shall only apply in the absence of mandatory legal regulations to the contrary.

KaVo cannot be held liable for defects and their consequences due to natural wear, improper cleaning or maintenance, non-compliance with operating, maintenance or connection instructions, calcination or corrosion, contaminated air or water supplies or chemical or electrical factors deemed abnormal or impermissible in accordance with factory specifications.

The warranty does not usually cover bulbs, glassware, rubber parts and the colour-fastness of plastics.

Defects or their consequences that can be attributed to interventions on or changes made to the product by the customer or a third party are excluded from the warranty. Claims from this warranty can only be asserted when the transfer form (copy) belonging to the product has been sent to KaVo, and the original can be presented by the operator or user.

1.5 Transportation and storage

1.5.1 Currently valid packaging regulations



Note

Only valid for the Federal Republic of Germany.

Properly dispose of and recycle the sales packaging, in accordance with the relevant packaging regulations, through waste management businesses or recycling companies that run a comprehensive return system. KaVo has licensed its sales packaging in accordance with this directive. Please conform with the regional, public waste-disposal system regulations.

1.5.2 Damage in transit

In Germany

If the packaging is visibly damaged on delivery, please proceed as follows:

1. The recipient of the package must record the loss or damage on the delivery receipt. The recipient and the representative of the shipping company must sign this delivery receipt.
2. Leave the product and packaging in the condition in which you received it.
3. Do not use the product.
4. Report the damage to the shipping company.
5. Report the damage to KaVo.
6. You must contact KaVo before returning a damaged product.
7. Send the signed delivery receipt to KaVo.

If the product is damaged but there was no discernable damage to the packaging upon delivery, proceed as follows:

1. Report the damage to the shipping company immediately and no later than 7 days after delivery.
2. Report the damage to KaVo.
3. Leave the product and packaging in the condition in which you received it.
4. Do not use a damaged product.



Note

Failure on the part of the recipient to comply with any of the above obligations will mean that the damage will be considered to have arisen following delivery (in accordance with ADSp. Art. 28).

Outside Germany



Note

KaVo is not liable for damage arising from transportation. Immediately inspect the delivery after receipt!

If the packaging is visibly damaged on delivery, please proceed as follows:

1. The recipient of the package must record the loss or damage on the delivery receipt. The recipient and the representative of the shipping company must sign this delivery receipt.
Without this evidence, the recipient will not be able to assert a claim for damages against the shipping company.
2. Leave the product and packaging in the condition in which you received it.
3. Do not use the product.

If the product is damaged but there was no discernable damage to the packaging upon delivery, proceed as follows:

1. Report the damage immediately or at least 7 days after the delivery to the delivery company.
2. Leave the product and packaging in the condition in which you received it.
3. Do not use a damaged product.



Note

Failure on the part of the recipient to comply with any of the above obligations will mean that the damage will be considered to have arisen following delivery (in accordance with CMR law, Chapter 5, Art. 30).






1.5.3 Information on the packaging: Storage and transportation



Note

Please keep the packaging in case you need to return the product for servicing or repair.

The symbols printed on the outside are for transportation and storage, and have the following meaning:

	Transport upright with the arrows pointing upwards
	Fragile - protect against impact.
	Protect from moisture.
	Permissible stacking load
	Temperature range

1 User instructions | 1.5 Transportation and storage

	Humidity
	Air pressure



2 Safety

2.1 Description of safety instructions

2.1.1 Warning symbol



	Warning symbol
---	----------------



2.1.2 Structure



	 DANGER
	<p>The introduction describes the type and source of the danger. This section illustrates the potential consequences of non-observance.</p> <ul style="list-style-type: none"> ▶ The optional step covers necessary measures for avoiding hazards.

2.1.3 Description of danger levels

Safety instructions with three hazard levels are used in this document for avoiding personal and property damage.

	 CAUTION
	<p>CAUTION indicates a hazardous situation that can cause damage to property, or mild or moderate physical harm.</p>

	 WARNING
	<p>WARNING indicates a hazardous situation that can cause death or serious injury.</p>

	 DANGER
	<p>DANGER indicates the maximum hazard level. indicates a hazardous situation that can directly cause death or serious injury.</p>

2.2 Purpose – Proper use



<p>Note Only the INTRAsurg 300plus with the IPX8 foot switch is permitted for use in operating rooms.</p>
--

The ambient temperature must lie between +10°C and +35°C.
The relative humidity must be between 30% and 75%.

This KaVo product is intended only for use in the field of dentistry, for surgery to expose and dissect oral tissue structures (such as the periodontal gap, gingiva, bone, the jaw, extractions and implantations). The product may not be used for a purpose for which it was not intended.

"Proper use" includes following all the instructions for use and ensuring that all inspections and service tasks are performed.

The overarching guidelines and/or national laws, national regulations and the rules of technology applicable to medical devices for start-up and use of the KaVo product for the intended purpose are to be applied and complied with.

The user must ensure that that the unit works properly and is in a satisfactory condition before each use.

During use, national legal regulations must be observed, in particular:

- the applicable health and safety regulations
- the applicable accident prevention regulations

Users have a duty to:

- Only use equipment that is operating correctly
- to protect himself, the patient and third parties from danger.
- to avoid contamination from the product

To guarantee constant readiness for use and maintenance of value of the KaVo product, the recommended annual servicing must be done. Yearly safety inspections are required.

KaVo recommends an annual service check. In this service check, the safety checks are performed according to *IEC 62353:2007* as well as calibration.

The safety check involves the following tests: Visually inspect the medical device and accessories; check the ratings of fuses that are accessible from outside; protective conductor check and leakage current measurements as per IEC 62353:2007; medical device function test with reference to accompanying documentation.

Authorized to repair and service the KaVo product:

- Technicians with appropriate product training from KaVo branches
- Technicians of authorized dealers specially trained by KaVo

In Germany, operators, equipment managers and users are obliged to operate their equipment in accordance with the MPG regulations.

The services encompass all the test tasks required in accordance with § 6 of the medical devices operator ordinance (Medizinprodukte-Betreiberverordnung, MPBetreibV).

After servicing, manipulation or repair of the device, the device must be tested according to IEC 62353:2007 (state of the art) before it is used again.

Check the parameters of the device settings after servicing. Enter the service parameters in the table on the last page.



Note

The product must be cleaned and serviced according to instructions if it is not to be used for a long period.



Note

Any waste which is generated must be recycled or disposed of in a manner which is safe both for people and for the environment. This must be done in strict compliance with all applicable national regulations. Questions on proper disposal of the KaVo product can be answered by the KaVo branch.



Note

A recycling pass can be downloaded from www.kavo.com.

2.2.1 Information on electromagnetic compatibility

Conformance level: The immunity test levels required in IEC 60601 are met.

	CAUTION
	<p>Damage from incorrect accessories</p> <p>The use of other accessories, transformers and lines than those indicated (with the exception of transformers and lines that KaVo sells as replacement parts for internal components) can increase transmission or reduce the electromagnetic immunity of the product.</p> <ul style="list-style-type: none"> ▶ Only use accessories recommended by KaVo.

	CAUTION
	<p>Interactions between devices</p> <p>Effect on proper use</p> <ul style="list-style-type: none"> ▶ Do not use the device directly next to or stacked with other devices ▶ If it is necessary to use it next to or stacked with other devices, check for proper use in this set up.



Note

Based on EN 60601-1-2 concerning the electromagnetic compatibility of electro-medical devices, we need to point out that:



- medical electrical devices are subject to special measures regarding electromagnetic compatibility and must be operated in accordance with KaVo assembly instructions.
- portable and mobile high-frequency communications devices can influence medical electronics.



Note

KaVo cannot guarantee that accessories, lines and transformers not delivered by KaVo will correspond with EMC requirements of EN 60601-1-2.

2.2.2 Risks from electromagnetic fields (implanted life-sustaining system)

	 DANGER
	<p>Risks from electromagnetic fields. The functions of implanted systems (such as pacemakers) can be influenced by electromagnetic fields.</p> <ul style="list-style-type: none"> ▶ Patients and users must inquire about implanted systems before using the motors and check the use. ▶ Evaluate the risks and benefits. ▶ Do not bring the motors close to the systems. ▶ Take suitable emergency precautions and immediately react to any changes in health. ▶ Symptoms such as elevated heart beat, irregular pulse and dizziness are signs of problems with a pacemaker.

2.2.3 Disposal of electronic and electrical devices



Note

According to EC directive 2002/96 concerning used electrical and electronic devices, this product is subject to the cited directive and must be disposed of accordingly within Europe.

Before disassembling / disposing of the product, it must be completely processed (disinfected, sterilized) according to the section "Preparation methods"

Additional information can be obtained from KaVo (www.kavo.com) or your dental supplier.

For final disposal, contact:

Germany







To return an electrical device, proceed as follows:


1. At the homepage www.enretec.de of enretec GmbH, you can download a form for a disposal request under the menu item `eom`, or you can use it as an online request.
2. Fill out the request with the corresponding information, and send it as an online request or by fax (+49(0)3304 3919 590) to enretec GmbH.
The following avenues are also available for questions and for initiating a disposal request:
Telephone: +49 (0) 3304 3919 500
E-mail: pickup@eomRECYCLING.com and
Post: enretec GmbH, eomRECYCLING Department
Kanalstraße 17
16727 Velten
3. Your **movable** device will be picked up in your practice, and your **permanently installed** unit will be picked up at the curb at your address on the agreed deadline. The owner or user of the device will bear the costs for disassembly, transportation and packaging.


International (EU)

For country-specific information on disposal, contact your dental supplier.

2.3 Safety instructions

	<p>⚠ CAUTION</p>
	<p>Electrical sparks in the product Explosion.</p> <ul style="list-style-type: none"> ▶ Do not use product in areas subject an explosion hazard.
	<p>⚠ CAUTION</p>
	<p>Damaged network cable/missing protective conductor. Electrical shock.</p> <ul style="list-style-type: none"> ▶ Check the network cable before use. The socket outlet must have a protective contact and meet the respective national guidelines. ▶ Always completely inset the power cable into the device's socket to connect to the mains.
	<p>⚠ CAUTION</p>
	<p>Damage due to liquids Faults on electric components.</p> <ul style="list-style-type: none"> ▶ Protect product openings from penetration of liquids. ▶ Remove liquids from the inside of the device.
	<p>⚠ CAUTION</p>
	<p>Unintended penetration of liquids. Electrical shock.</p> <ul style="list-style-type: none"> ▶ Do not immerse the product in a tub-like container. ▶ Check and make sure that the coolant containers and lines do not leak before each use. If liquid is found on the device, do not touch it and disconnect the power cable from the mains. Completely dry the surface of the unit before re-inserting the power cable.
	<p>⚠ CAUTION</p>
	<p>Damage to the instrument hoses from stickers. Instrument hoses can explode.</p> <ul style="list-style-type: none"> ▶ Do not affix stickers or adhesive tape.
	<p>⚠ CAUTION</p>
	<p>Rotating parts while the pump is operating Injury</p> <ul style="list-style-type: none"> ▶ Do not stick anything in the pump. Turn off the device when the pump is open.


	⚠ CAUTION
	<p>Risks from electromagnetic fields. The functions of implanted systems (such as pacemakers) can be influenced by electromagnetic fields.</p> <ul style="list-style-type: none">▶ Ask patients before treatment!

	⚠ CAUTION
	<p>Electrostatic charge. Destruction of electronic components.</p> <ul style="list-style-type: none">▶ Before contacting the hose coupling, touch the bottle holder with your hand to release static electricity.


KaVo recommends that only **original KaVo parts®** be used for operating and repairs since their safety, operation and specific suitability have been tested in extensive tests.

3 Product description

3.1 INTRAsurg 300/300 plus

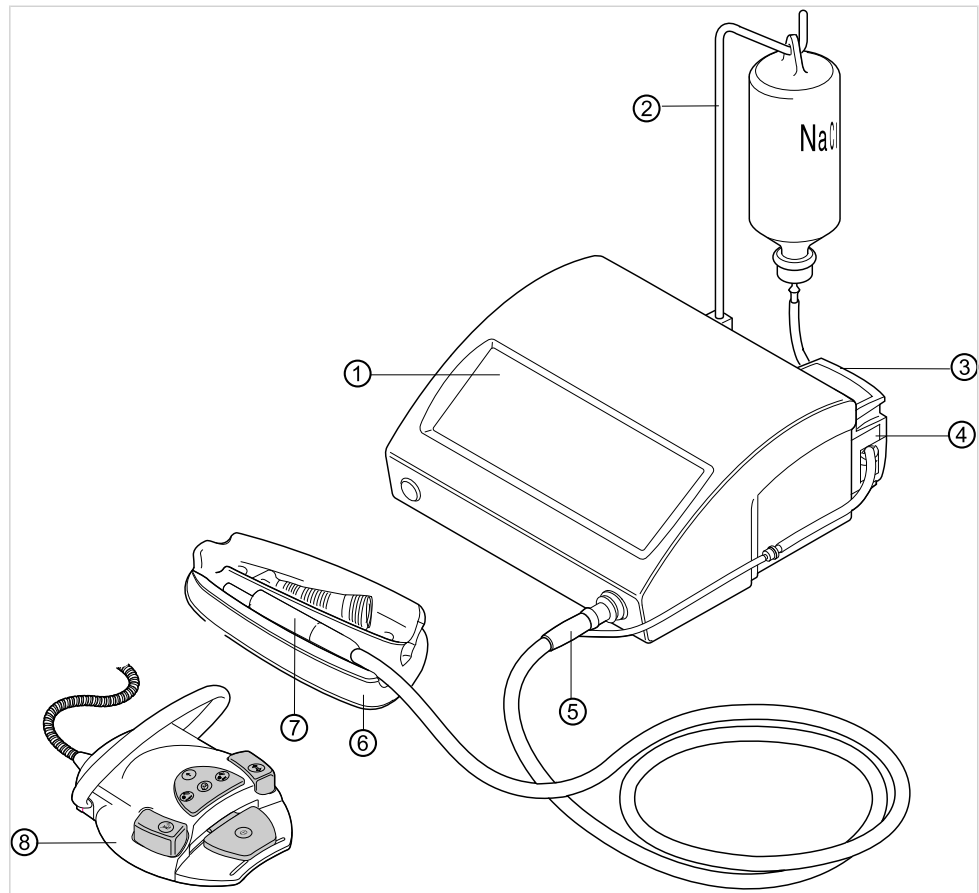

- | | |
|--------------------------------|--------------------------------|
| ① Controls | ⑤ Surgical hose coupling |
| ② Bottle holder | ⑥ Motor holder |
| ③ Unlocking the pump mechanics | ⑦ Surgery motor |
| ④ Hose pump | ⑧ Multifunctional foot control |

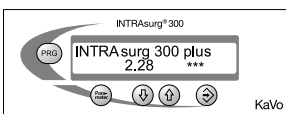


The KaVo INTRAsurg 300 plus has an additional lighting function and recognition function for the INTRA CL light, hand and contra-angle handpiece.



Note

Additional information on surgery and implantology can be obtained at: www.ka-vo.com

3.2 Controls

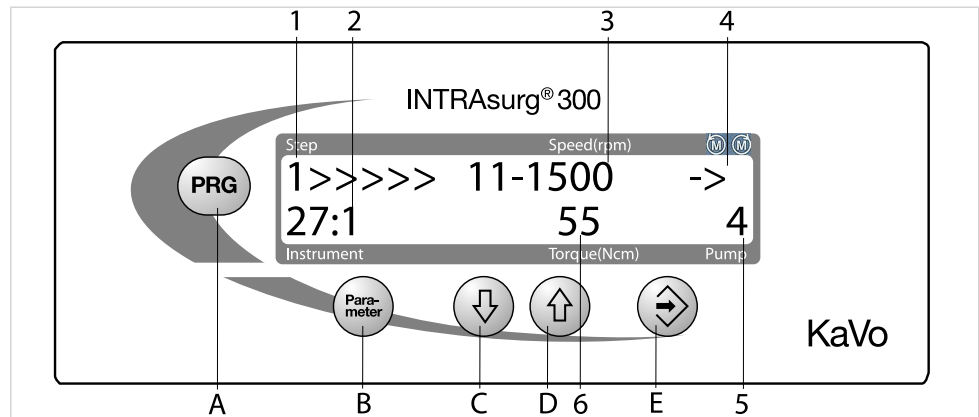

- A Program button
- B Parameter button
- C Minus button
- D Plus button
- E Enter button

- ① Step (1... 6)
- ② Transfer ratio of the handpiece and contra-angle handpiece
- ③ Speed (rpm)
- ④ Direction of rotation (< - / - >)
- ⑤ Coolant pump setting (1... 6)
- ⑥ Maximum torque (Ncm)




Open the program. Use the PRG button to open the steps; repress the button to switch to free use.



Open setting mode. And Use the parameter button to select changeable parameters. Press the parameter button until the asterisk (*) flashes at the desired position.

Changeable parameters:

- max. torque
- max. speed
- coolant pump setting
- Transfer ratio of the handpiece and contra-angle handpiece
- Direction of motor rotation

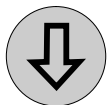

Reduce * value.



Increase * value.



Enter button:

Save set of values/leave setting mode


INTRAsurg 300 Plus: Request maximum torque.



Note

The motor can only be started after you leave the setting mode by pressing the ENTER button. This saves all the changes. Only after this is done is the motor released to start.

3.3 Multifunctional foot control

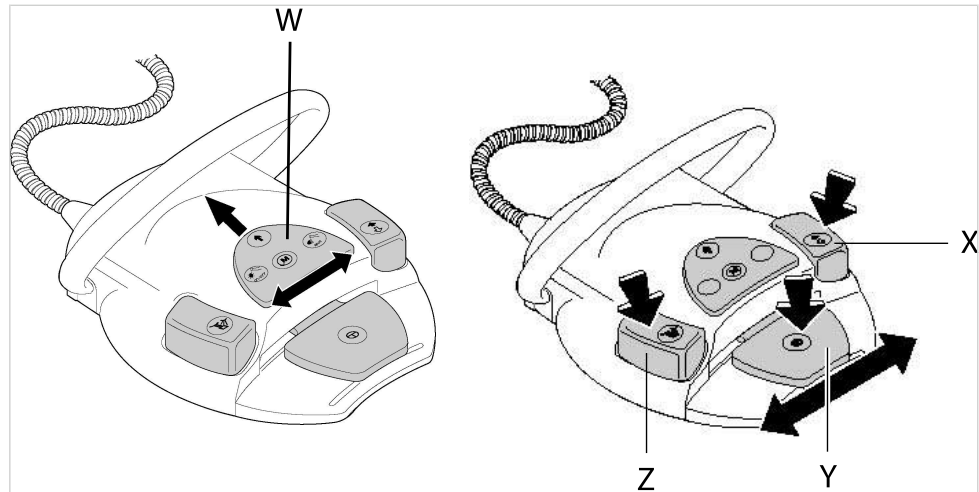

- (W) Shift pedal
- (X) Right button

- (Y) Hinged pedal
- (Z) Left button

Shift pedal W:

Move forward for the default counter-clockwise rotation of the surgical motor
 Only with the INTRAsurg 300 Plus: The light can be switched on/off by moving to the left, and the spot light can be turned on by moving to the right.

Right button X:

Procedural step button for selections in the program and for advancing the individual steps. Press briefly = proceed forward; press for a long time = proceed backwards. Double-click = open selection. An acoustic signal confirms proceeding forward/backward when the motor starts.

Hinged pedal Y:

Footswitch to start the motor and control its speed, and to start withdrawing the coolant without the motor running and then and then regulate the amount conveyed. The process is started by pressing the pedal down. The intensity is decreased by pressing the pedal to the left, and increased by pressing it to the right.




Note

Always check the input values before each use.

Left button Z:

Coolant button to turn the coolant on/off and select rinsing without the motor.

3.4 Rating plate for 300/300 plus

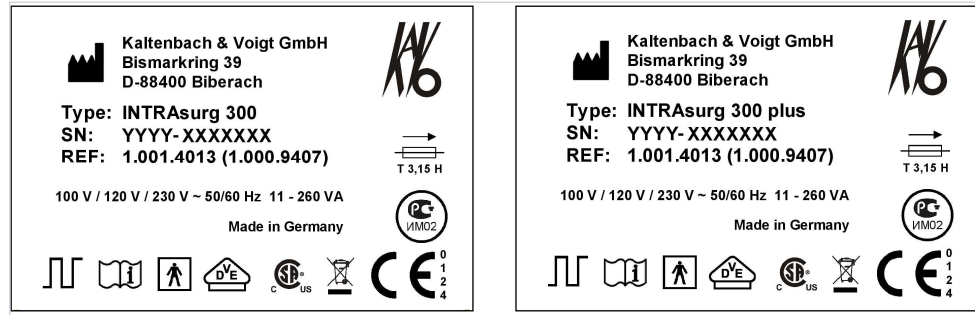

	Installation side and securing the mains input
	CE mark
	VDE mark
	CSA mark
	Classification, type BF
	Follow instructions for use
	Mode: Continuous operation with intermittent load
	Supply voltage
SN:	YYYY = Year manufactured XXXXXXX = Serial number
REF:	Material number
Type:	Device type
	For disposal information, see use in accordance with intended purpose
	GOST R certification


3.5 Technical data for the 300/300 plus

Width	291 mm
Depth	306 mm
Height	126 mm
Weight of the device	6.85 kg
Weight of foot control	2.35 kg
Weight of motor	approx. 170 g

Input voltage	100 V~/120 V~/230 V~
Input frequency	50/60 Hz
Power consumption	Max. 260 VA
True power	Approximately 7 Watt
Speed	300 – 40000 rpm
Pump delivery rate	25 – 150 ml/min
delivery pressure	1.5 – 2.2 bar
Multifunctional foot control	Protection class IPX1
Only KaVo INTRAsurg 300 Plus: Multifunctional foot control, identification: blue FBE stirrup, rating plate. Permissible in surgery rooms.	Protection class IPX8
Standard-length hose for the multifunctional foot control	1.90 m
Special-length hose for the multifunctional foot control Only available with the INTRAsurg 300plus.	3.5 m
Standard motor hose length Mat. no. 1.001.2651	2.00 m
Special length motor hose Mat. no. 1.004.6825	3.0 m
Mode continuous operation duty type	30 sec. operation/ 9 min.sec pause
The torque precision of the KaVo angle piece 27:1 (TRA C09 + C3) ranges from 20 to 30 Ncm at 20 to 50 rpm. Greater deviations are possible with other angle pieces.	± 3 Ncm
Transportation and storage conditions	
Ambient temperature	-20°C to +70°C
Relative humidity	5% to 95%
Air pressure	700 hPa to 1060 hPa

4 First use

4.1 Electrical connection

	⚠ CAUTION
	<p>Damaged network cable/missing protective conductor. Electrical shock.</p> <ul style="list-style-type: none">▶ Check the network cable before use. The socket outlet must have a protective contact and meet the respective national guidelines.▶ Always completely inset the power cable into the device's socket to connect to the mains.


Both lines L and N are individually fused with T 3.15 H. The network frequency can be 50 or 60 Hz.

Requirement

Make sure that the supply voltage is correct.

- ▶ First plug the inlet connector for non-heating apparatuses into the device, and then plug the other end of the line into a socket.

See also: 4.3 Adjust the mains input voltage, Page 19

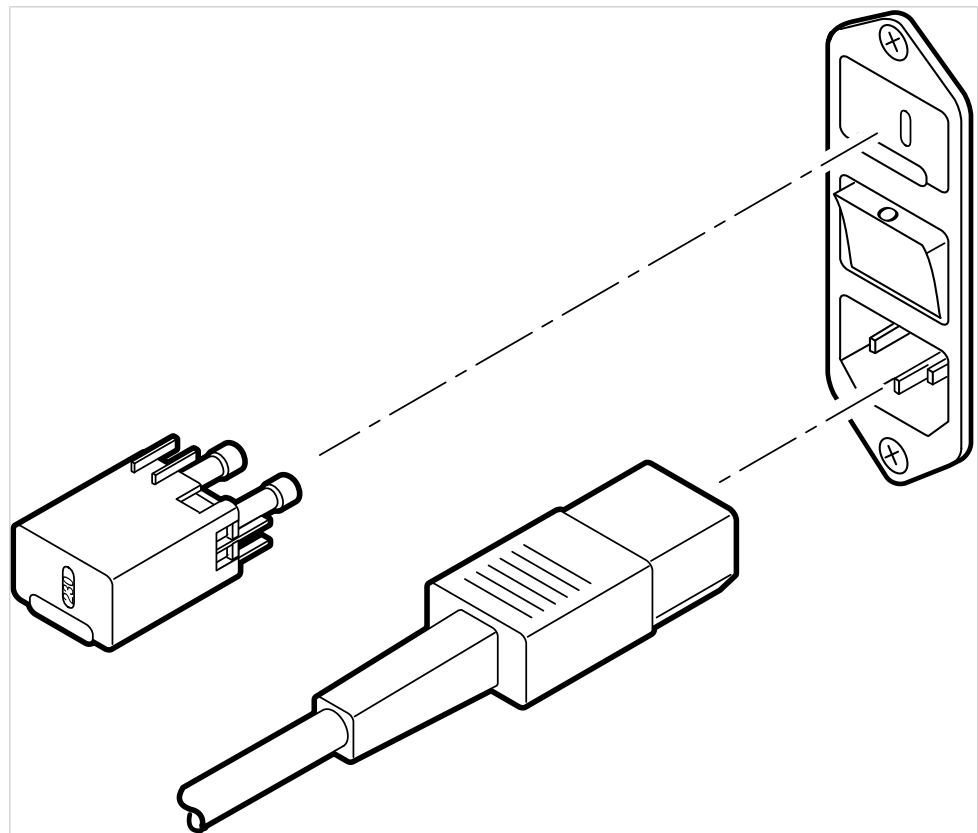

4.2 Unpacking



Note

Do not kink the motor hose since it may damage it.



Note

Save the box with all the packaging materials so that the device can be safely shipped (for example for a yearly service check). Always send both hose boxes so that the device can be safely packaged to prevent damage.


Prevent damage by properly removing from the package.

- ▶ Open box.
- ▶ Remove the hose box. Remove the multifunctional foot control and additional equipment.



Note

When removing, note that the multifunctional foot control is tightly connected to the device.

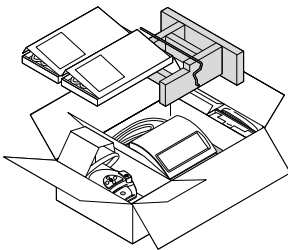

- ▶ Pull out the INTRAsurg 300/300 plus vertically upward and place on a flat surface.

4.3 Adjust the mains input voltage



⚠ CAUTION

Voltage-conducting parts

Electrical shock.

- ▶ Unplug the power cable from the power supply.


A country-specific mains input voltage (a rating play) is set at the plant. Check the setting and change it if necessary to the existing mains input voltage. The voltage selector ① can be set to 100V, 120V or 230V. The existing setting can be read at the cutout ④.

Remove the drawer ③.

Remove the fuse holder ② from the fuse drawer ③, and turn it so that the required voltage value appears in the cutout ④.


4 First use | 4.3 Adjust the mains input voltage

Shove the fuse holder ② back in the drawer ③. Insert the drawer ③ in the device.
A code protects the drawer from being incorrectly inserted.

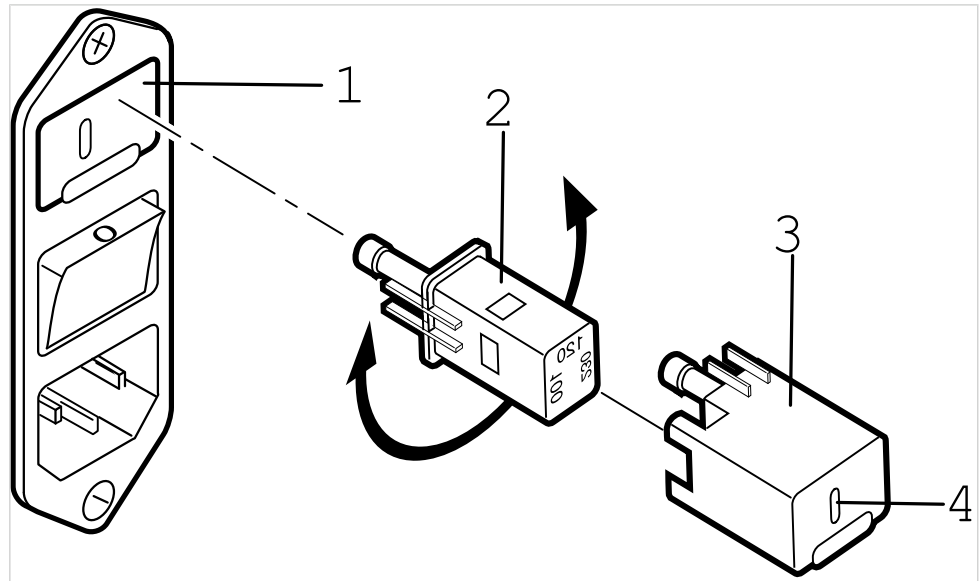

4.4 Mount the bottle holder and use the correct coolant container

- ▶ Plug-in the bottle holder ① **Mat. no. 0.761.1872** and rotate it until the bottle holder faces away from the device and always secure with screw ② **Mat. no. 1.005.8527**.

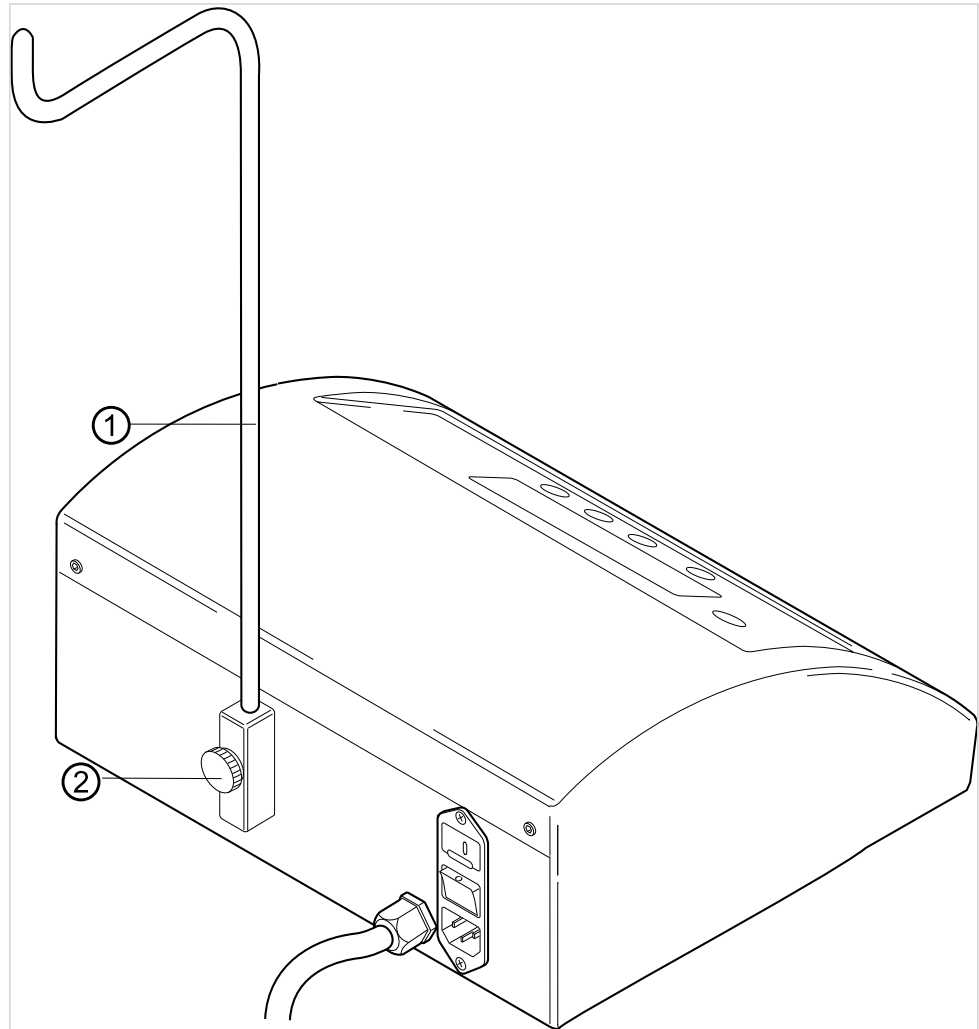

⚠ CAUTION

Tipping danger from excessively heavy coolant reservoir
Malfunctions

- ▶ Do not use coolant reservoirs that can hold more than 1 litre!
- ▶ Ensure that the device is standing securely.

In the case of glass bottles, insert an aeration needle. Otherwise, the supplied amount decreases with prolonged removal of coolant.



⚠ WARNING

Electric shock due to coolant ingress

Injuries


- ▶ Rotate the bottle holder away from the device to ensure that no leaking coolant drips on the device or the electrical connection.

4.5 Connect the surgery motor and coolant reservoir



Note

The delivered parts are not sterile. Sterilise them before the first patient treatment. All parts conducting liquids must be sterilised.

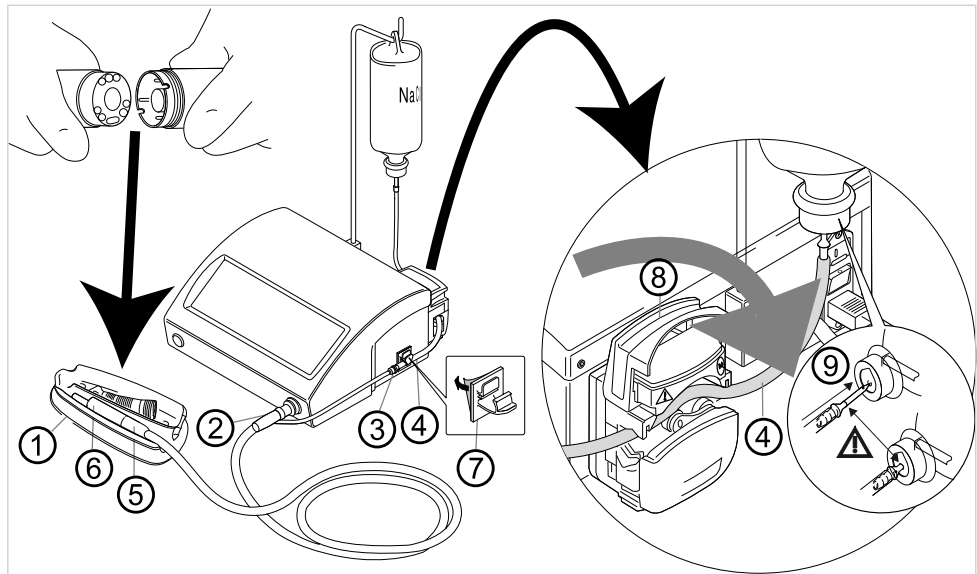


- ▶ Plug in the surgery motor ⑥ at these surgery hose coupling and tighten with a union nut ⑤. Follow the separate instructions for use for the motor.
- ▶ Place the surgery motor on the motor holder ①.




Note


The kink-protection part on the coolant hose inlet may not be twisted. Twisting the kink protection can cause defects in the hose and lead to a restriction of coolant flow and causes the coolant to drip afterwards.

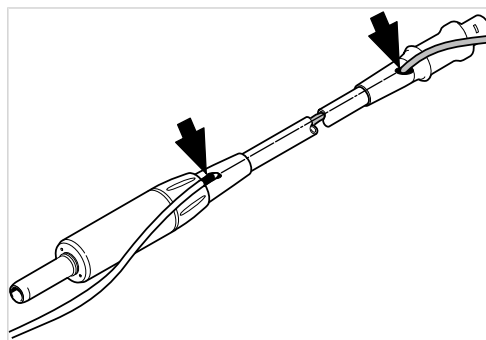



ESD symbol

	<p>⚠ CAUTION</p>
	<p>Electrostatic charge. Destruction of electronic components.</p> <ul style="list-style-type: none"> ▶ Before contacting the hose coupling, touch the bottle holder with your hand to release static electricity.

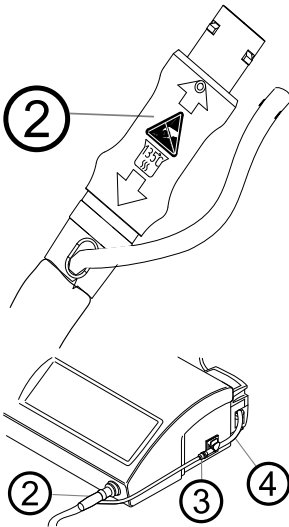


- ▶ Insert the hose coupling ② into the device connection until it audibly locks into place and the marking ring on the plug is covered.
- ▶ Shove the coolant hose ④ onto the plug-in nipple ③ by twisting slightly.
- ▶ Open the pump mechanics by lifting the lock ⑤.




Note

Insert the pump hose into the pump so that it is not clamped or pinched by the lock. Run all hoses relaxed without tension. Shove the hoses firmly on the plug-in nipple to ensure a tight connection.

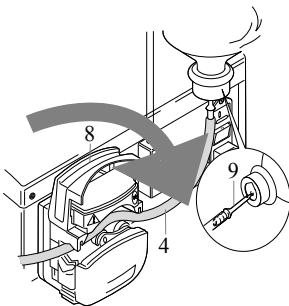

- ▶ Close the lock ⑧ in the direction of the arrow.
- ▶ Insert the insertion needle ⑨ into the coolant reservoir, and hang the coolant reservoir on the bottle holder. Check the seal and seat of the insertion tip. Prevent fluid from exiting the device.


	⚠ CAUTION
	<p>Tipping danger from excessively heavy coolant reservoir Malfunctions</p> <ul style="list-style-type: none"> ▶ Do not use coolant reservoirs that can hold more than 1 litre! ▶ Ensure that the device is standing securely.

Insert a sterile ventilation needle into glass bottles. The flow will otherwise be reduced if coolant is withdrawn for a long period. Check the level of the coolant reservoir several times after waiting a while. The system must be leak-free. No liquid may escape from the coolant hose or the side of the connections.



Note

The saline solution corrodes all the metal parts. Immediately remove residue.

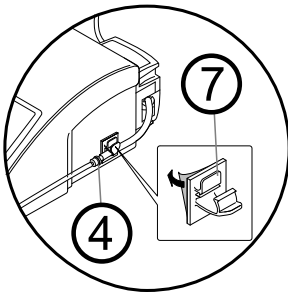

- ▶ Fasten the self-adhesive cable clip ⑦ **Mat. no. 1.001.3348** to the INTRAsurg 300/300 plus (for the location, see the illustration).


- ▶ Insert the coolant hose ④ into the cable clip and close.



Note

Always monitor the level of the coolant reservoir. The device itself does not monitor the level. Use transparent coolant reservoirs.


If the product is used on potentially infectious patients, use disposable products for supplying coolant and thorough infection protection. KaVo recommends using a sterile disposable protective hose over the motor and hose (approx. 1.2 m) and a transparent disposable protective film over the device. The control buttons and display must always be identifiable for operation.

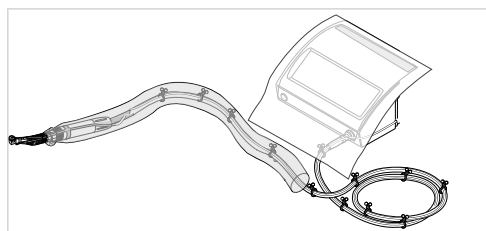


5 Operation

5.1 Operation in general

5.1.1 Start up

	⚠ CAUTION
	<p>Product damage from incorrectly adjusting the tension. Destruction of electrical components.</p> <ul style="list-style-type: none"> ▶ Adjust the voltage input for the product to the existing mains voltage.

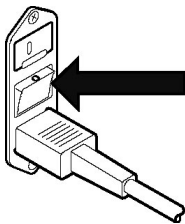

- ▶ Switch on the device.

An acoustic signal confirms that the device is turned on.
The device runs a self test.
The firmware version and note "Please check parameters" is briefly shown on the display.

- ▶ Note the displays.


If the device is not monitored, we recommend turning it off for safety reasons and to save energy.



Note

The last saved values are displayed.

5.2 Surgical motors

	⚠ CAUTION
	<p>Incorrectly set parameters Damage from incorrect input values</p> <ul style="list-style-type: none"> ▶ Check all input values before use.




Note

Following instructions for use, service instructions and installation instructions in the motor, handpiece and contra-angle handpiece packaging.

5.2.1 INTRA S 550 surgical motor

- ▶ Plug straight or contra-angle handpiece onto the motor coupling.
- ▶ Check secure seat before each treatment.

- ▶ Attach coolant hose.

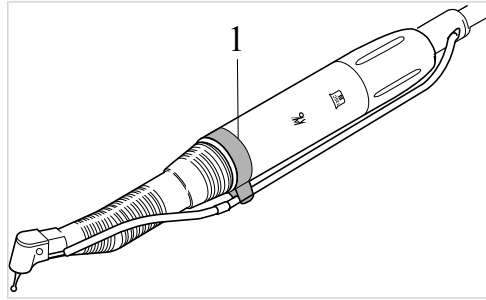


5.2.2 INTRA LUX SL 550 surgical motor



Note

The INTRAsurg 300 plus has a recognition function for INTRA CL light, hand and contra-angle handpieces. Handpiece and contra-angle handpiece-specific parameters are automatically transferred. Check all input values before use.

Mount the handpiece or angle piece

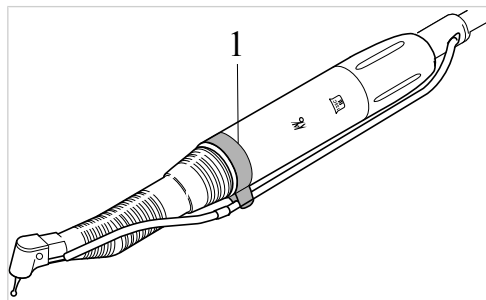



⚠ CAUTION


Heat transferred from high-pressure.
Burns from touching the high-pressure lamp.

- ▶ Only operate the motor with the retention ring ① mounted.


- ▶ Place the **INTRA CL** light, hand and contra-angle handpiece on the surgery motor and turn it until the catch audibly locks in place.

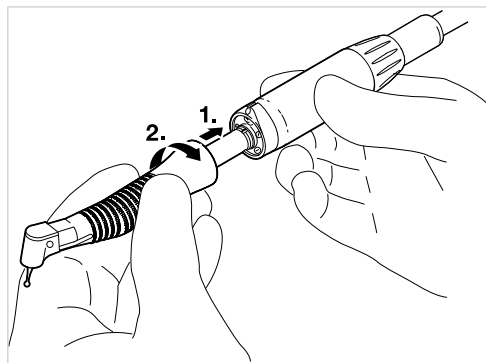

	⚠ CAUTION
	<p>Damage from changing the handpiece and angle piece during operation. Wear to the catch on the handpiece and angle piece and the motor Unbalanced motor axis.</p> <ul style="list-style-type: none"> ▶ Change the handpiece and angle piece only when the motor not running.

	⚠ CAUTION
	<p>Damage from operating with open chucks The handpiece or angle piece locks and rotates on its own axis.</p> <ul style="list-style-type: none"> ▶ Only operate the handpiece or angle piece when the chuck is locked.

5.2.3 Remove the handpiece or angle piece

	⚠ CAUTION
	<p>Damage from changing the handpiece and angle piece during operation. Wear to the catch on the handpiece and angle piece and the motor Unbalanced motor axis.</p> <ul style="list-style-type: none"> ▶ Change the handpiece and angle piece only when the motor not running.




Note

Prevent to coolant hose from dripping on the motor. Coolant may not penetrate the motor.

- ▶ Remove the coolant hose from the handpiece or angle piece.
- ▶ Twist the handpiece or contra-angle handpiece slightly while pulling it off.

See also: 5.1.1 Start up, Page 25

5.2.4 Start the motor

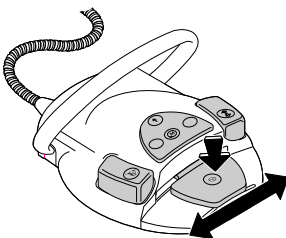

- ▶ Start the motor. Press the slide pedal and change the speed by moving it the side. Left stop: minimum speed; right stop: maximum preselected speed. When the slide pedal is no longer actuated, the selected speed is set on the multifunctional foot control. Repress the slide pedal to retrieve the last selected speed.


See also:

5.3.3 Preset the speed range, Page 30
5.4.4 Set parameters, Page 36



Note

When the hinged pedal is pressed in an unchanged position (to the side), the last selected position with an assigned motor speed is retrieved. The motor immediately runs at the assigned speed. After the program step changes, check the speed setting, and move the hinged pedal to the required position before pressing.

	⚠ CAUTION
	<p>Excessive speed. Injuries.</p> <ul style="list-style-type: none"> ▶ In applications where the speed is critical, move the hinged pedal to the left position (minimum speed) before the motor starts and move up to the desired speed. ▶ Set the speed to a safe maximum for the application.

5.3 Free use

5.3.1 Parameters menu


The following values can be changed in the parameter menu:

- Transfer ratio of the handpiece and contra-angle handpiece
- max. speed
- Direction of motor rotation
- Coolant pump setting
- max. torque




- ▶ Press the parameter button.

A star ① flashes next to the alterable parameter.

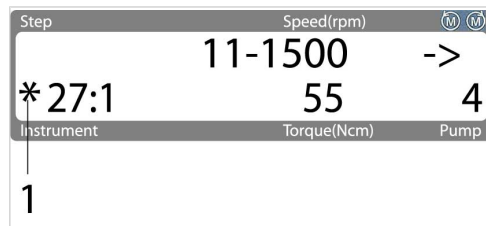

- ▶ Select the desired parameter by pressing the parameter button several times.

An asterisk * flashes next to the alterable parameter.

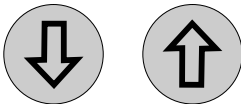

- ▶ Set the parameters to the desired value with the plus or minus key.



- ▶ Press the Enter key.

The set values are saved.



Note

Always save the selected setting with enter button when leaving the menu. Otherwise, the motor will not start (safety circuit).

Examples of saving values of handpieces and contra-angle handpieces



Note

Default values for the respective transfer ratios can be selected with the Plus or Minus buttons.

Indication	Transfer ratio	Combined unit Reducing shank	Head	Speed range rpm	Maximum torque Ncm
Grinding, cutting	1:2	3555 1:2	--	600-80.000	2,7
Grinding, Drilling, cutting	1:1	3,610 N 1:1	--	300-40.000	5,5
Grinding, Drilling, Thread, Implant	12:1	3624 N 4:1	C 3 3:1	25-3.300	40
Grinding, Drilling, Thread, Implant	16:1	3624 N 4:1	67 RIC 4:1	18-2.500	45
Grinding, Drilling, Thread, Implant	27:1	C 09 9:1	C 3 3:1	11-1.500	55
Grinding, Drilling, Thread, Implant	36:1	C 09 9:1	67 RIC 4:1	9-1.100	55




Note

The listed indications are only examples. To prevent unnecessary risk, observe the guideline speeds given by the manufacturer of the rotating instruments. KaVo recommends setting a speed range of 11-25/min or 11-35/min for screwing tasks.

5.3.2 Set the torque limitation




- ▶ Press the parameter button until the asterisk flashes at the torque display ③.

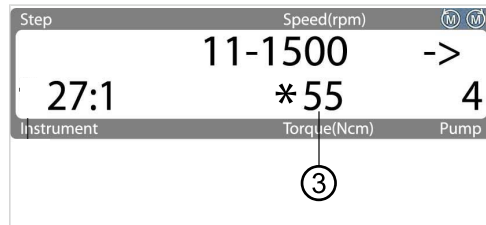

Note

The torque values ③ are only for KaVo handpieces and angle pieces that operate properly.

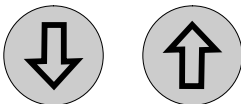

The maximum value can be changed in the parameter menu by using the Plus and Minus buttons.



Note

Once 90% of these set torque has been reached, a signal sounds. When the set torque is exceeded during treatment, the motor stops. The foot pedal has to be pressed again to start the motor.

Checking the handpieces and contra-angle handpieces:

Check the set transfer ratio, and let the motor and handpiece or contra-angle handpiece run for 20 seconds at maximum speed. The torque ③ is an indicator for the condition of the instrument.

Display of 0 = status is good

Display < 0.5 = improve maintenance

Display > 0.5 perform service

5.3.3 Preset the speed range





Note

Only the maximum value ① can be preset. Use the multifunctional foot control to select the speed within the predetermined range by simultaneously pressing the slide pedal and moving it to the side.



- ▶ Press the parameter button until the asterisk flashes at the maximum speed ①.

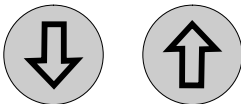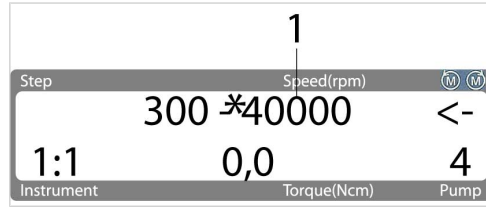



The maximum value can be changed in the parameter menu by using the Plus and Minus buttons.

	⚠ CAUTION
	<p>Excessive speed. Injuries.</p> <ul style="list-style-type: none"> ▶ In applications where the speed is critical, move the hinged pedal to the left position (minimum speed) before the motor starts and move up to the desired speed. ▶ Set the speed to a safe maximum for the application.

5.3.4 Set the coolant flow



- ▶ Press the parameter button until the asterisk flashes at the coolant flow ⑤.

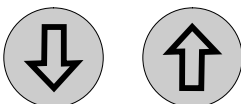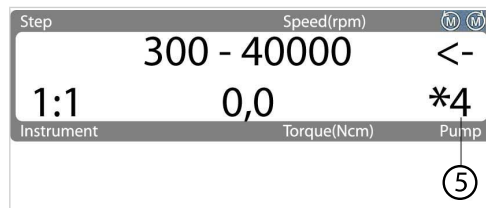


The flow ⑤ can be changed in the parameter menu by using the Plus and Minus buttons.

1: Pump at the minimum flow for approximately 25 mL/min

6: Pump at the maximum flow for approximately 150 mL/min

The flow is determined by the size of the outlet. When the outlet is small, back-pressure arises in the hose that can cause the hose connections to separate. Dripping coolant indicates that this may occur.

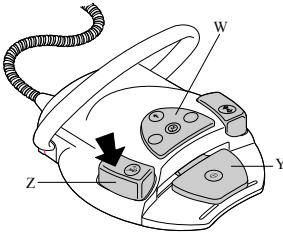



Note

The pump can only be turned on and off with the left button on the multi-functional for control. Display when the pump is turned off ⑤: 0

To save the pump at OFF (0), press the button on the multifunctional foot control and save it with the enter button. The previously set flow is retained. If you can turn it on again by pressing the left button on the multi-functional foot control and save it with the enter button.

5.3.5 Selecting the handpiece and contra-angle handpiece

- ▶ Press the parameter button until the asterisk flashes at the transfer ratio ③.

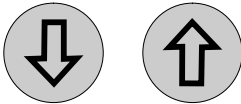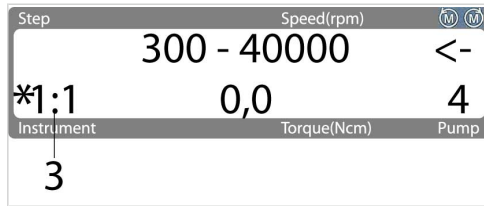

Use the Plus and Minus buttons to scroll through the specified table of transfer ratios ③ of the instruments in the parameter menu. Always check the selection in the display.



Note



The maximum permissible speed and torque of the selected handpieces and contra-angle handpieces are displayed.

5.3.6 Switching the direction of motor rotation



Note

The direction of motor rotation can only be switched with the multi-functional for control. Counterclockwise rotation <- ④ remains active until the next switch. If counterclockwise rotation is not saved by pressing the ENTER button, the motor rotates clockwise -> when the device is turned on. Wait for the acoustic anti-clockwise signal.

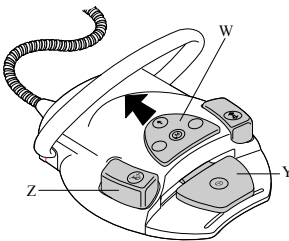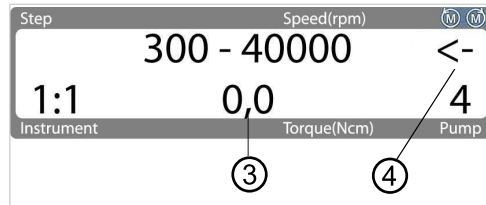


- Move the shift pedal in the direction of the arrow.

With counterclockwise rotation <- the rotational direction display ④ flashes, and two acoustic signals sound. When the motor is started, three acoustic signals sound. In the PRG menu, counterclockwise rotation can be saved by pressing the enter button.




Note

No display of the torque value ③ in counterclockwise rotation.

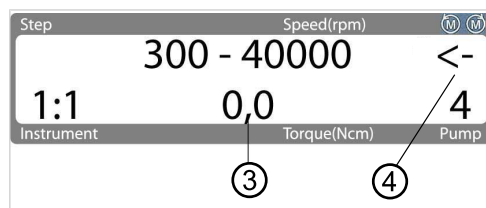


5.3.7 Rinsing function

This function can be selected from each step, but it cannot be saved.

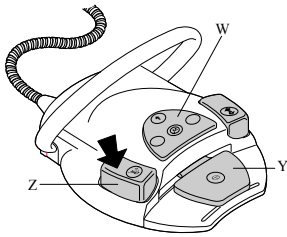


- ▶ Press the left button for 3 seconds to select the rinsing function.

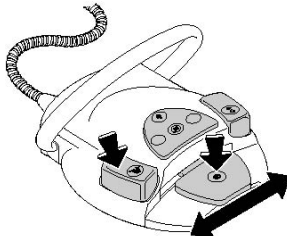

- ▶ Start the drawing of the coolant by pressing the slide pedal. The flow can be changed in six steps by moving the slide pedal to the side.



- ▶ The rinsing function is shown in the display.

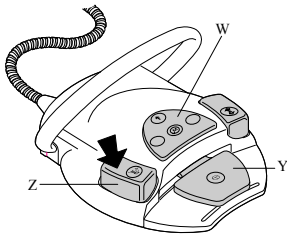

- ▶ Press the left button briefly to stop the rinsing function. The device returns to the previously selected function.

5.3.8 Surgical motor in free use



Note

This function is only supported by the **INTRAsurg 300 plus** in conjunction with the **INTRA CL** light, handpiece and contra-angle handpieces.


The recognition function of the **INTRA CL** light, handpiece and contra-angle handpieces enable parameters for the free use of handpiece and contra-angle handpieces to be recognised and adopted.

5.3.9 Request maximum torque



Note

This function is only supported by the **INTRAsurg 300 plus**.

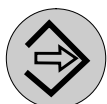

- ▶ Hold down the enter button after the motor stops.

The maximum torque ① of the last motor operation is displayed in Ncm.

Step	Speed(rpm)	M	M
	11 - 1500	->	
27:1	max 24	4	
Instrument	Torque(Ncm)	Pump	
	1		



Note

The torques are overwritten when the motor restarts under a minimum load of 5 Ncm.


5.4 Program

5.4.1 Start the program

1

Step	Speed(rpm)	M	M
1>>>>>	11 - 1500	->	
27:1	55	4	
Instrument	Torque(Ncm)	Pump	

The product has a basic program with six programmable steps. The current step is shown in the display ①.

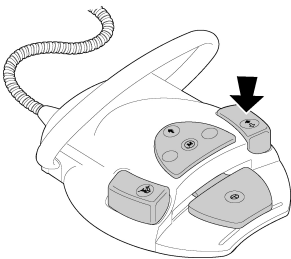

- ▶ Press the right button briefly.

or




- ▶ Press the program button once.

An acoustic signal confirms the program selection.




- ▶ To switch to free use from the program, press the program button again briefly.

5.4.2 Switch forward step

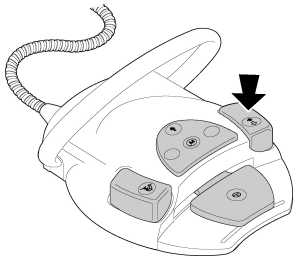

- ▶ Press the right button briefly.


or

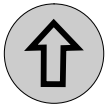

- ▶ Press the plus button.


An acoustic signal confirms the entry.

5.4.3 Go backward a step

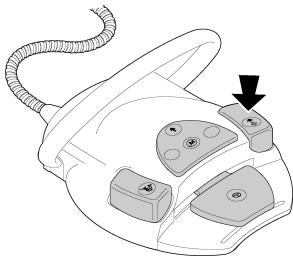

- ▶ Press the right button for a while.


or

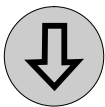

- ▶ Press the minus button.

An acoustic signal confirms the entry.

5.4.4 Set parameters

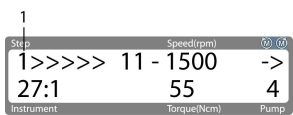

- ▶ Select the step ① to be changed.



- ▶ Switch to the parameter menu by pressing the parameter button.

See also: 5.3 Free use, Page 28


The values shown for a step are default values that you can use to start work immediately. They can all be changed in the parameter menu and can hence be adapted to the your individual approach.

Changed values can be saved in the parameter menu and are then available for the next use.

The following values can be changed in the parameter menu:

- Transfer ratio of the handpiece and contra-angle handpiece
- max. speed
- Direction of motor rotation
- Coolant pump setting
- max. torque

See also: 5.3 Free use, Page 28


	⚠ CAUTION
<p>Incorrect use. danger of injury.</p> <ul style="list-style-type: none"> ▶ Always check the values before use. 	

5.4.5 Limit the number of steps

Six steps are saved upon first start up and after setting up the delivery status.

Example: Limiting the steps

The number of steps should be limited to four.

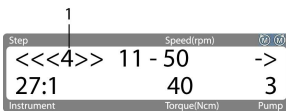

▶ Select step 4 ①.



▶ Press the program button for at least three seconds.


An acoustic signal confirms the entry.



Note

Now only the first four steps can be selected.
When you continue to step 4, the message **END** appears.
If you continue, you go to free use.

Remove restriction

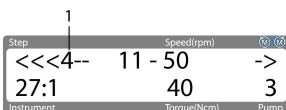

▶ Select the last step ① (4 in this example).



- ▶ Press the program button for at least three seconds.

An acoustic signal confirms the entry.

5.4.6 Setting up delivery status



Note

All entries and changes made are deleted with this function.



Note

If the device was initialised, everything saved previously is deleted. Therefore, for simplicity, everything that was saved should be noted down for easy recovery.



- ▶ Turn device on while pressing the programme button.

After initialization, the device is ready for use in its condition at the time of delivery.

Condition at the time of delivery: Programme

Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Direction of motor rotation
1	27:1	11-1500	55	4	->
2	27:1	11-50	40	4	->
3	27:1	11-50	40	4	->
4	27:1	11-50	40	4	->
5	27:1	11-50	55	0	-<
6	27:1	11-50	10	0	->

Condition at the time of delivery: Free application

Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Direction of motor rotation
	27:1	11-1500	55	4	->

5.4.7 Surgical motor in the program

Requirement

The mounted **INTRA CL** light, hand and contra-angle handpiece must be selected on the display.



If the mounted handpiece or contra-angle handpiece is not selected, the device asks for the corresponding input.

if the motor is turned on without the corresponding entry, three acoustic signals sound, and the read-in instrument flashes on the device display with question marks, for example: ?27:1? .




Note

The motor will not start (safety circuit).
To start the motor, either the selected handpiece or contra-angle handpiece must be mounted, or the mounted handpiece or contra-angle handpiece must be saved using the ENTER button.


Save the light, hand and contra-angle handpieces (without the mounted and peace and contra-angle handpiece)



- ▶ Press the parameter button until the asterisk flashes at the transfer ratio.

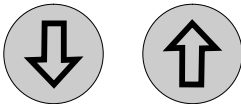


- ▶ Set the correct parameters with the Plus or Minus button.

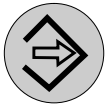

- ▶ Save the set values with the enter button.

The data are now saved in this step.

5.4.8 Request maximum torque

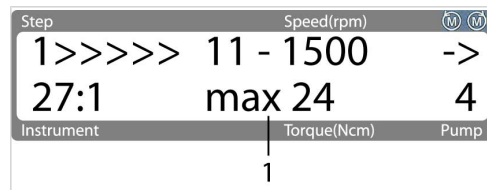


The maximum torque is saved when you advance a step.
When the motor is started in the next step, the highest torque is always saved.



- ▶ Press the enter button after closing the program.

In each step ①, the max. torque ② is shown in the display.

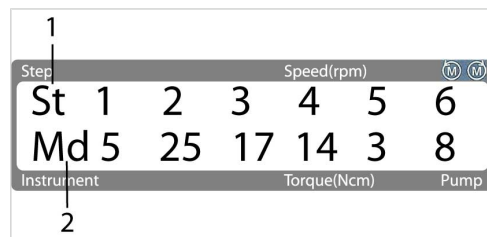



- ▶ To leave the display, press the enter button.

All saved torques are deleted.



Note

If you leave the program by pressing the plus button or the right button on the multifunctional foot control, all torques are deleted. If you accidentally go too far and skip the torque display, the values can be reactivated and displayed by pressing the right button on the multifunctional foot control.

5.5 Functions exclusively with the 300 plus program version

Additional functions of the INTRAsurg 300 plus with an INTRA LUX SL550 motor:


- Light function.
- IPX8 multifunction foot control (blue stirrup, rating plate).
- Recognition of the INTRA CL light, hand and angle piece. The transfer ratio and associated values are shown on the display.
- The user can perform additional calibration with the INTRAsurg Calibration 1.002.3569. (27:1 angle piece is required)
- Display of the maximum torques. The maximum torque is saved in each step. The max. torque can be read on the display during surgery or after the program ends (only with clockwise rotation).
- Read-out table for the torque values at the end of applications.

5.5.1 Recognition function




Note


To prevent misinformation, the contact's provided for recognition always need to be regularly cleaned on the motor and instrument pusher. The contacts always must be shiny. Wipe off the contacts with a cloth soaked in disinfectant and rub dry.

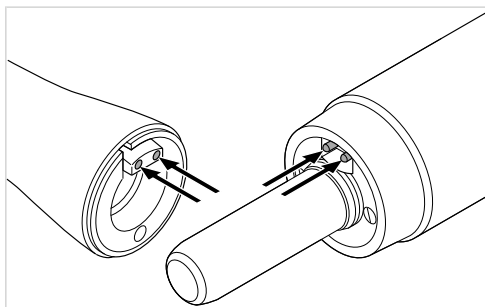


- ▶ Place the INTRA CL light, hand and angle handpiece on the INTRA LUX SL 550 motor and turn it until the catch ① audibly locks in place.

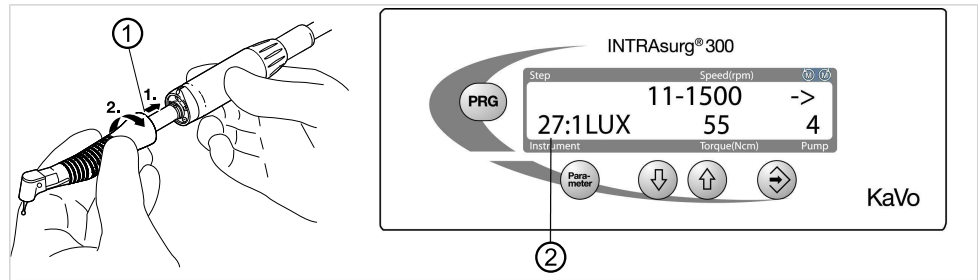

The hand and angle piece is recognised, and the corresponding saved parameters are shown on the device display ②.

	<p>⚠ CAUTION</p>
	<p>Incorrect input values. Danger of injury.</p> <ul style="list-style-type: none"> ▶ The input values must be checked before each use.

These parameters are immediately activated in free use.

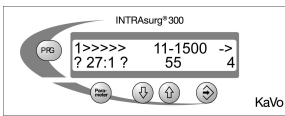


The recognised handpiece or angle piece must be selected in the program. If it is not selected, the device asks for the corresponding input. Three acoustic warning signals sound, and the read-in instrument flashes on the device display with a question mark, for example: ?27:1?



Note

For safety reasons, the motor does not start (safety circuit) since the input value on the display always has priority.

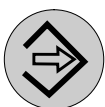

- ▶ Check the parameters. Values displayed by the device can be saved by pressing the ENTER button.


The data of the mounted INTRA CL light, hand and angle piece are only saved in this step. The motor starts with these released parameters.



Note


This must be done for the saved parameters to be reliably maintained and continuously followed.

5.5.2 Light functions

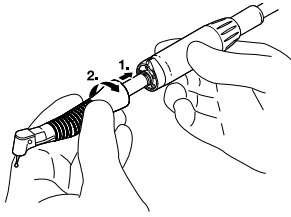

Place the INTRA CL light, hand and contra-angle handpiece on the INTRA LUX SL 550 motor and turn it until the catch audibly locks in place. The INTRA CL light, hand and contra-angle handpiece is recognised and the corresponding parameters appear on the display (this function is only possible with the INTRA CL light, hand and contra-angle handpiece).


Light on/off without motor and pump running (spot)

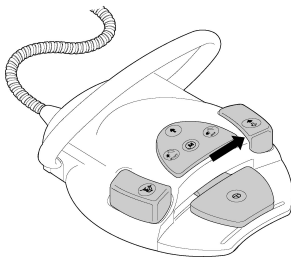

- ▶ Move the shift pedal to the right.

The light only shines when actuated (spotlight function).

Default light when the motor starts

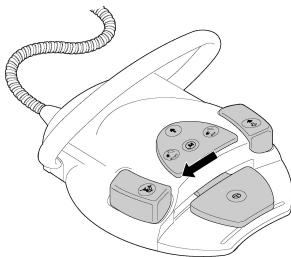

- ▶ Move the shift pedal to the left.

The light is turned on or off. When the light is turned on, the following appears in the display: LUX (this function is only possible with INTRA CL light, handpieces and contra-angle handpieces).



Note

If there is no light function, check if the light is turned on and if the INTRA CL light, hand and contra-angle handpiece is properly mounted.


Afterglow time

The light shines for 3 sec. after the motor stops running.



Note

The motor can only operate with retention ring ① **Mat. no. 1.001.5029.**

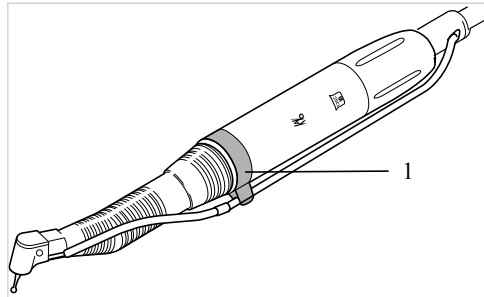

Note

Only the light function is activated and a 1:1 hand and angle piece is displayed with hand and angle pieces with a metal snap-in button. The torque is not displayed.

5.5.3 Request maximum torque (free use)



- ▶ Hold down the enter button after the motor stops.

The maximum torque ① of the last motor operation is displayed in Ncm.

Step	Speed(rpm)	M	M
	11 - 1500	->	
27:1	max 24	4	
Instrument	Torque(Ncm)	Pump	
	1		



Note

The torques are overwritten when the motor restarts under a minimum load of 5 Ncm.

5.5.4 Request maximum torque (program)

See also: 5.4.8 Request maximum torque, Page 39


5.5.5 Recommended programming when setting a number of implants one after the other

Mode of operation:

Perform all tasks for drilling an implant in Free use mode. Screw in the respective implant in PRG steps (max. 6). At the end of the task, the torque screw-in value can be requested and manually documented for evaluating primary stability (in the form of the following table, for example).

Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
1	27:1	11 – 50	40	0	->
2	27:1	11 – 50	40	0	->
3	27:1	11 – 50	40	0	->
4	27:1	11 – 50	40	0	->
5	27:1	11 – 50	40	0	->
6	27:1	11 – 50	40	0	->

5.5.6 INTRAsurg calibration

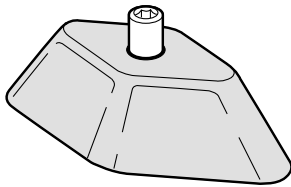

The influence of the efficiency of the instruments is taken into account when testing the torque display.

Store the INTRAsurg calibration device at 15°C to 30°C to provide precise calibration. The calibration value refers to the temperature of the INTRAsurg calibration of 23°C. If this temperature is not maintained, the value needs to be corrected according to table. For example: Calibration at 20°C (read 15.2) results in a corrected input value (+0.3) of 15.5.


Temperature	15°C	16°C	17°C	18°C	19°C	20°C	21°C	22°C	23°C	24°C	25°C	26°C	27°C	28°C	29°C	30°C
correction	+0.8	+0.7	+0.6	+0.5	+0.4	+0.3	+0.2	+0.1	0	-0.1	-0.2	-0.3	-0.4	-0.5	-0.6	-0.7

A sterile instrument can also be used for calibration. To prevent contamination, use a sterile catch pin.



Note

The INTRAsurg Calibration system can be disinfected.

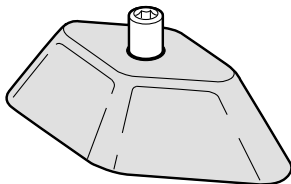


- ▶ Pre-position cleaned and set up angle piece 27:1 and INTRAsurg Calibration Mat. no. 1.002.3569.



- ▶ Press the parameter button for 3 seconds.



The following message appears on the display:

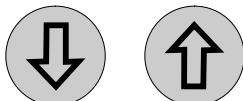


- ▶ Read the calibration affixed to the bottom of the INTRAsurg Calibration (such as 15.2), correct the value according to the temperature, and enter with the Plus and Minus keys.



- ▶ Confirm with the enter button.

- ▶ Insert the catch pin into the INTRAsurg Calibration.
- ▶ Place the INTRAsurg Calibration on a fixed base and hold.

- ▶ Start the motor. The motor always rotates to the right.

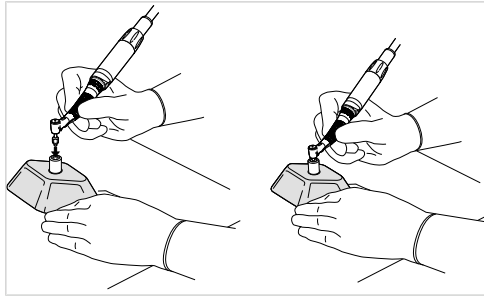

Note

Do not apply any pressure.
The motor may not experience any additional load.

The process can be monitored on the display.



After calibration, the following appears on the display:



If calibration is unsuccessful, the following appears on the display:



Note

The output values are retained in case of interruption or error.

The calibration automatically ends. The INTRAsurg 300 plus is ready with the new setting.

6 Preparation methods DIN EN ISO 17664



Note

Before preparations, follow the instructions for use of the surgery motor, handpieces and contra-angle handpieces, and the scaler SONICflex and special tools (bone).

The instructions for cleaning and sterilisation have been validated by the manufacturer. Any deviation from the provided instructions by the person preparing the device should be checked to see if it is effective, and potential negative consequences should be evaluated.

Sufficient safety that prevents microorganisms that have entered the instruments from becoming a possible cause of infection is only provided by carefully processing the outside and inner surfaces directly after use. The handpieces and angle pieces must be packaged and sterilised and sterile when used.

If potentially infectious liquids and particles can contact the products, it is recommendable to cover and protect these areas with sterile disposable products.

Frequent preparation with proper use does not significantly affect the unit and motor. The product life normally ends due to wear and damage from use. Observe all manufacturer instructions.

Treat KaVo gently when cleaning them. Remove deposits such as saline coolant residue, blood and saliva with a disposable cloth immediately after use. Only use fully functional products for treating patients.

Observe all current national hygiene regulations. See also www.rki.de (infection protection)



Note

A German language CD on hygienic measures in medical practices is obtainable from: MEDIEN & Forum, Am Brunnen 8, 29229 CELLE – Tel. D 05141/370188/89. E-mail: hammer@medien-information-bildung.de. Additional information can be obtained at: www.rki.de

Site of use

Remove moisture-sensitive parts close to the disinfection site. The treatment surfaces of the area should be hygienically smooth surfaces. Protect difficult to clean surfaces and objects from contamination by covering them.

Preparation for decontamination



- Wear suitable protective clothing.
- Wear eye protection. No eating, drinking or smoking.
- Disinfect hands properly.
- Turn off the device and unplug it.
- Prepare and follow a hygiene plan.
- Only use disinfectants suitable for this group of devices.
- Prepare the parts suitable for sterilisation: Care, rinsing and packaging.
- Use qualified trained personnel for preparing dental equipment and instruments. Train the personnel according to these instructions for use Plan and provide a schedule.

6.1 Cleaning

Always clean directly after using the medical device to prevent contaminants from hardening. Saline solution contains sodium chloride which corrodes metal. Quickly remove all residue to prevent surface damage.

6.1.1 Machine cleaning

Not applicable for INTRAsurg device parts.

	 CAUTION
	<p>Damage due to liquids Faults in electric components.</p> <ul style="list-style-type: none"> ▶ Do not machine-clean the INTRAsurg device.

6.1.2 Manual cleaning

- ▶ Wipe off the supply unit (all visible surfaces), motor hose, motor holder, bottle container, foot control surfaces and the connecting lines with a damp disposable cloth. Do not wipe the cloth on edges and press so that liquid will not penetrate. Remove all residue and dry the device.



Note


After uses that run the risk of infection, rinse the coolant lines of the handpieces and contra-angle handpieces with a disinfectant (approx. 20 ml) with a disposable syringe. Avoid damage. Rinse the coolant lines with demineralised water, then drain and dry them.

6.1.3 Rinse the hoses

Preparation for product decontamination:


Directly after each use, rinse all the hoses to prevent malfunctions from the crystallisation of the saline solution. All of the hoses must be completely emptied. Then sterilise them and observe the hygienic requirements for storage and maintaining sterility.

- ▶ Remove the pump hose needle from the coolant reservoir. Immerse the pump hose needle in a clean container with demineralised water (50-60 ml). Place the spray hose from the instrument in an empty container or sink.

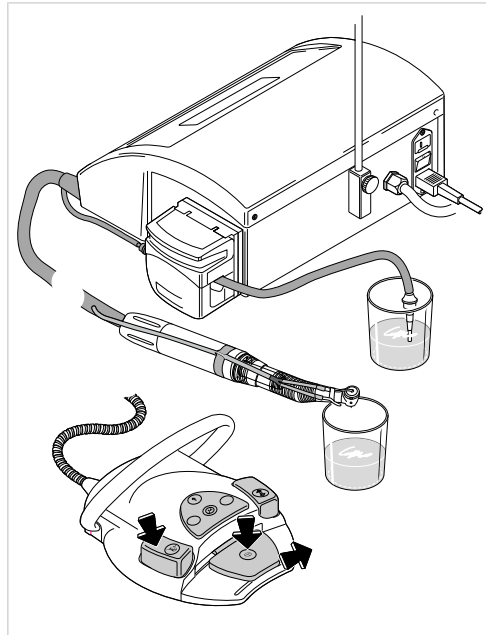

Note

Rinse until all the hoses are drained. Do not immerse the water outlet into the liquid. Observe the necessary measures for infection protection.

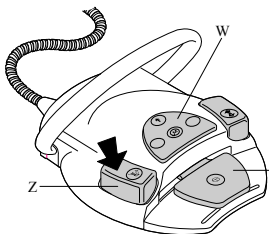


- ▶ Press left button 3 sec. The rinsing function is selected.


See also: 5.3.7 Rinsing function, Page 33



- ▶ The rinsing function is shown in the display.

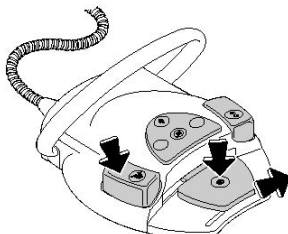

- ▶ Hold down slide pedal and slide it to the right until the stop. The hose pump runs at the max. rate.
- ▶ Ensure that the water jet is in the right direction and that it is the right amount.

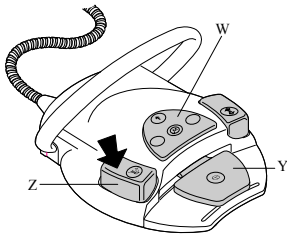

- ▶ Press the left button briefly. You exit the rinsing menu.

6.2 Disinfection



Note


The motor and motor hose cannot be thermidisinfected. After treating each patient, disinfect the surfaces close to the patient that are contaminated by contact or aerosols. Wipe disinfection should be used in all cases.

6.2.1 Automated disinfection

Not applicable for INTRAsurg device parts.

6.2.2 Manual disinfection

- ▶ Pull the insertion needle ③ out of the coolant reservoir ②.
- ▶ Open the lock ① and remove the hose ④.

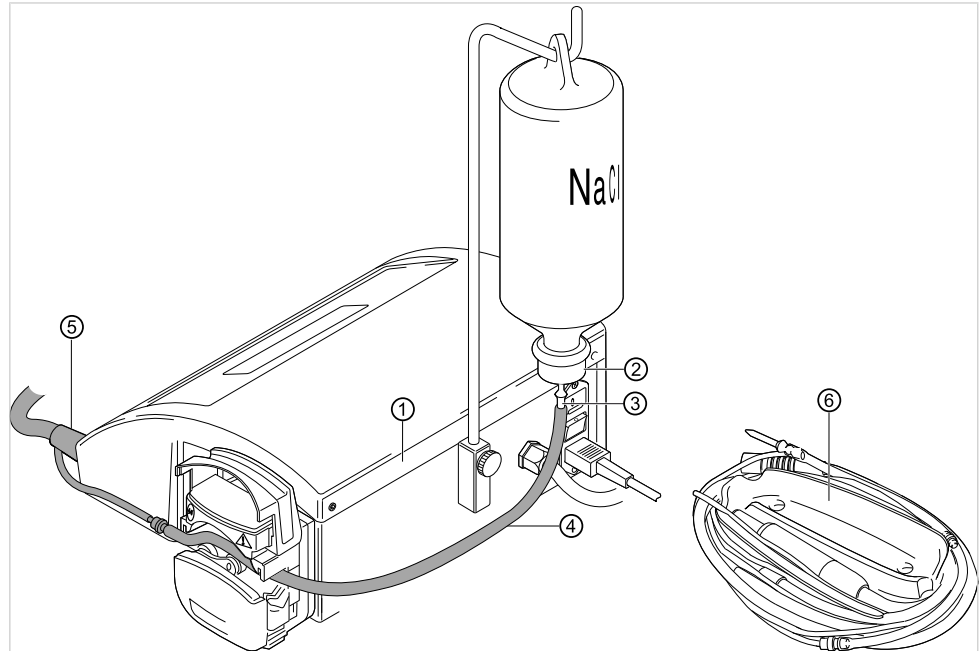

- ▶ Remove all hoses marked gray including tips and surgical motor on the unit.
- ▶ Wipe-disinfect surfaces with a soft, disposable cloth and permissible disinfectants. Make sure that the device is completely wet. Let the disinfectant act for the prescribed time. Dry the surface.

Permissible disinfectants (the uses correspond to the existing manufacturer's instructions and national guidelines. observe the safety data sheets.) KaVo recommends the following products based on material compatibility. The microbiological efficacy must be ensured by the disinfectant manufacturer.

- Incidin Liquid (Ecolab)
- FD 322 Dürr
- Schülke & Mayr Microcide
- Isopropanol 70%




6.2.3 Drying


- ▶ After disinfection and sterilisation: Let the surfaces and plug-in parts for the hoses and cables (contacts) completely dry in the room before they are used again.


6.2.4 Service, inspection and testing after preparation

- ▶ Test of cleanliness and intactness, care, repair: check the hose connections for ease of plugging and secure attachment. Check the function of device settings and motor operation. Check hose pump for sufficient coolant. Check secure attachment of hose connections. Check control commands on foot control. Check recognition function of INTRAsurg 300 plus. Comply with hygiene requirements (sterility) during the check. If there are any sites of breakage or clearly recognised changes on the surface, the parts must be checked by the Service.

6.3 Sterilisation in a steriliser in compliance with DIN EN 13060

	<p style="background-color: yellow; padding: 2px;">⚠ CAUTION</p> <p>Product damage due to improper sterilisation. Damage to the sterile product.</p> <ul style="list-style-type: none"> ▶ No hot air sterilisation, no chemical cold sterilisation, do not sterilise with ethylene oxide!
	<p style="background-color: yellow; padding: 2px;">⚠ CAUTION</p> <p>Product damage due to improper sterilisation. Damage to the sterile product.</p> <ul style="list-style-type: none"> ▶ An initial vacuum can shorten the product life.
	<p style="background-color: yellow; padding: 2px;">⚠ CAUTION</p> <p>Moisture Non-sterility</p> <ul style="list-style-type: none"> ▶ Ensure dryness. Autoclaves with a subsequent vacuum ensure dryness. In addition, drying can be accelerated with a 10 minute drying phase with the autoclave door open.

	⚠ CAUTION
	<p>Product damage Contact corrosion</p> <ul style="list-style-type: none"> ▶ Remove the sterilised item from the autoclave immediately after sterilising and drying.

	⚠ CAUTION
	<p>Crystallising saline solution Malfunctions</p> <ul style="list-style-type: none"> ▶ Rinse all hoses and drain them completely before they are sterilised.



Note

The user is responsible for observing the regulations and conditions for sterility. After each patient, dispose of the coolant reservoir, rinse the hoses, and completely drain and sterilise them.



Note

When treating patients who may have an acute, critical infectious disease, be sure to observe the hygienic measures cited in applicable publications and reports. If possible, use suitable disposable products to avoid the transmission of critical pathogens. These protect the user, the patient and all participants in the surgery. A microbiologically appropriate coolant liquid must be used. Only use prescribed isotonic saline solution E NaCl 0.9 (identified as an infusion solution; follow the instructions on the package insert) for cooling and spraying wounds. All materials that are contaminated from the dental and medical fields must suitably processed and sufficiently identified after cleaning and sterilisation (RKI recommendation, www.rki.de).

See also: 6.1.3 Rinse the hoses, Page 47




The KaVo medical devices released for sterilisation have a maximum temperature resistance of 138 °C.

6.3.1 Packaging

The sterilisation bag must be large enough for the items to be sterilised so that the bag is not stretched. The quality and use of the packaging of the items to be sterilised must satisfy the applicable standards and be appropriate for sterilising.

- ▶ Wind all hoses previously cleaned and disinfected that are marked grey including tips and the surgical motor on the motor holder ⑤.
Seal the motor holder, handpieces and contra-angle handpieces in sterile bags or place them on a surgical tray.

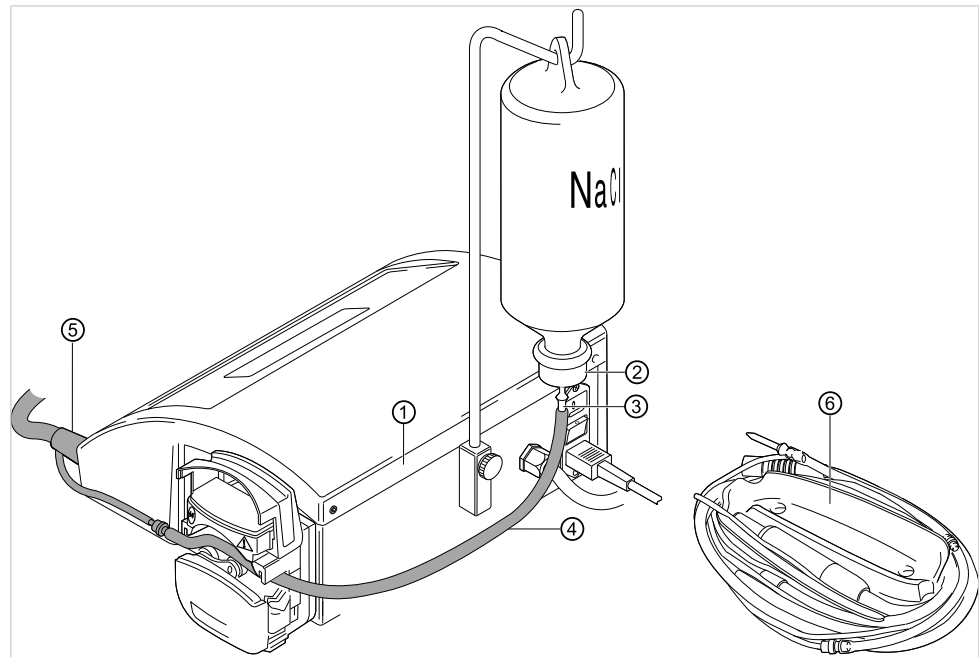

6.3.2 Sterilisation

- ▶ Sterilise in saturated steam at $135\pm 1^{\circ}\text{C}$, gravity, for at least 4 minutes.



Note

Let the sterilised materials cool and try to the ambient temperature before using them again.

6.3.3 Storage

Observe all necessary measures for hygiene when storing sterile goods. Store protected from dust, release with identification on the packaging. Monitor storage length.

7 Troubleshooting



Note

If your problem cannot be rectified after reading the following troubleshooting documentation, please contact a KaVo-trained Service Technician.

In case of a problem, a fault number is displayed with: "**warning xx**" or "**remark xx**"
Fault numbers 1 to 67 are possible.



Note

Confirm all other error messages with the ENTER button. If the error message does not disappear or the error is reported again, please notify customer service.
Please always write down what you want or the service you want performed when you send in the device. The settings must be rechecked after the device is returned.

Malfunction	Cause	Remedy
Nothing works.	Fuse blown.	▶ Check and if necessary replace the fuse.
	The unit must be switched off.	▶ Turn on the mains switch on the back of the unit.
	Both ends of the power line are not plugged in.	▶ Plug in the mains power line:
	The wrong handpiece or angle piece is mounted or set on the device.	▶ Adjust the right handpiece or angle piece on the device.
	Unknown.	▶ Turn the machine off and on.
No coolant in the instrument.	Coolant flow not preselected.	▶ Set the coolant flow. See also: 5.3.4 Set the coolant flow, Page 31
	The pump hose is sticky or crusty.	▶ Exchange the pump hose. See also: 7.1 Exchange the pump hose, Page 56
The coolant in the instrument is insufficient.	The spray nozzles are crusty or soiled.	▶ Clean the spray nozzle with the nozzle needle.
	The pump hose is sticky or crusty.	▶ Exchange the pump hose. See also: 7.1 Exchange the pump hose, Page 56
	Glass bottle with saline solution is not ventilated.	▶ Ventilate the glass bottle.
Coolant is dripping	Clogged coolant line or restricted flow	▶ Check if the coolant hose and connecting nipple are blocked.

Malfunction	Cause	Remedy
	The coolant pump is set too high.	<ul style="list-style-type: none"> ▶ Adjust the coolant flow. <p>See also: 5.3.4 Set the coolant flow, Page 31 7.1 Exchange the pump hose, Page 56</p>
	The pump hose is too long.	<ul style="list-style-type: none"> ▶ To avoid dripping, keep the pump hose length after the pump as short as possible.
Leaky hose pump housing.	The pump hose is worn.	<ul style="list-style-type: none"> ▶ Exchange the pump hose. <p>See also: 7.1 Exchange the pump hose, Page 56</p>
Leaky surgical hose.	Worn or defective coolant hose.	<ul style="list-style-type: none"> ▶ Exchange coolant hose. <p>See also: 7.2 Exchange coolant hose, Page 57</p> <ul style="list-style-type: none"> ▶ Shorten the hose by approx. 6 mm.
The motor shuts off during treatment without an error message.	The wrong handpiece or angle piece is mounted or set on the device.	<ul style="list-style-type: none"> ▶ Adjust the right handpiece or angle piece on the device.
The motor makes a grinding noise.	The motor is not correctly plugged on or screwed on.	<ul style="list-style-type: none"> ▶ Firmly insert the motor hose into the housing.
		<ul style="list-style-type: none"> ▶ Firmly screw on the motor hose to the motor.
		<ul style="list-style-type: none"> ▶ Turn the machine off and on.
		<ul style="list-style-type: none"> ▶ Check if all the connections and couplings are firmly seated.
The motor runs out of true.	The motor is not correctly plugged on or screwed on.	<ul style="list-style-type: none"> ▶ Firmly insert the motor hose into the housing.
		<ul style="list-style-type: none"> ▶ Firmly screw on the motor hose to the motor.
		<ul style="list-style-type: none"> ▶ Turn the machine off and on.
		<ul style="list-style-type: none"> ▶ Check if all the connections and couplings are firmly seated.
No light on the handpiece and angle piece.	The light is not turned on.	<ul style="list-style-type: none"> ▶ Turn on the light.

Malfunction	Cause	Remedy
	The hand and angle piece is im- properly attached.	▶ Mount the handpiece and angle piece until the catch audibly locks.
	The lamp is faulty.	▶ Replace the lamp Mat. no. 1.002.2928. See also: Instructions for use for the motor.
	The catch is dirty.	▶ Clean the catch (gold contacts). See also: Instructions for use for the motor and instrument.
	No suitable light, hand and angle piece	▶ Use a suitable light, hand and angle piece.
Display: >> remark 65 -motor is overloaded	Motor overload warning	▶ Stop treatment as soon as pos- sible and let the motor cool off for a few minutes.
Display: >> warning 50 -motor is overloaded	Motor is overloaded, the maximum permissible motor output was ex- ceeded.	▶ The motor is blocked for 10 mi- nutes; let it cool.
Display: >> remark 66 -motor not plugged	The motor hose is not correctly in- serted.	▶ Properly insert the motor hose.
	The motor is not correctly connec- ted to the motor hose.	▶ Properly plug in the motor.
Display: >> remark 67 -motor not plugged.	When the motor starts, the discon- nection of the motor from the hose is detected.	▶ Plug the hose onto the motor and screw it tight. This display disappears after eight seconds. The motor can be restarted after this time.
Display: >> warning 36 -missing foot control	Defective footswitch or defective connection to the footswitch, possi- bly by unintentional penetration of moisture.	▶ Dry the footswitch. Turn the machine off and on
Display: error 3	Internal system error	▶ Exchange the pump hose. Con- tact service.
Display: error 4	Internal system error	▶ Exchange the pump hose. Con- tact service.
Display: error 29	The motor hose is damaged.	▶ Check the motor hose and ex- change it if necessary
Display: warning 49	The power consumption of the sur- gical motor was too high.	▶ Press the ENTER button and continue. Exchange the motor hose. Exchange surgical motor.

7.1 Exchange the pump hose



Note

The parts delivered for exchange are not sterile. Rise and sterilize them before the first treatment.


See also: 6.3 Sterilisation in a steriliser in compliance with DIN EN 13060, Page 50



Note

If the pump hose becomes brown after extended use, it must be exchanged.

- ▶ Make sure that all lines containing saline solution are completely empty.
- ▶ To drain the lines, remove the insertion needle ③ from the saline solution reservoir ② and activate the pump by using the multifunctional foot control until the hoses are free of saline solution. Remove the container and store it so that it does not drain.

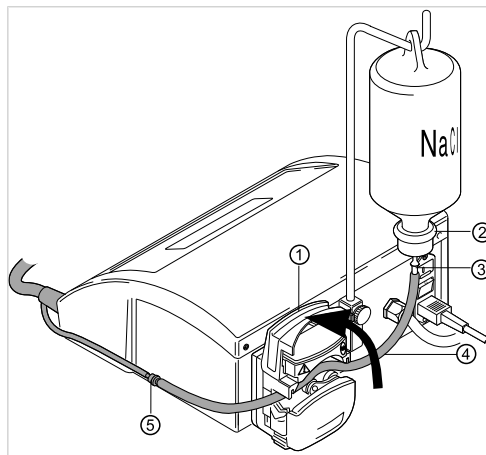

- ▶ Turn off the INTRAsurg 300/300 plus.




⚠ CAUTION

Running, open hose pump.
Injury hazard

- ▶ Turn off the device before opening the hose pump.

- ▶ Lift the lock and remove the pump hose.
Insert the new pump hose **Mat. no. 1.001.3347**.

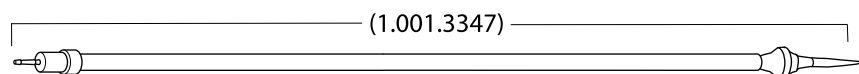



Note

Insert the pump hose into the pump so that it is not clamped or pinched by the lock. Run all hoses relaxed without tension. Shove the hoses firmly on the plug-in nipple to ensure a tight connection.



Note

To avoid dripping, keep the pump hose length after the pump as short as possible.


- ▶ Carefully turn the lock ① downwards and turn it until it locks into place.

The pump hose length after the pump determines the dynamic pressure in the coolant hose after the motor or pump stops (drop formation).

Pump hose ① short: less accumulated volume in the coolant hose (less potential dripping)

Pump hose ② long: large accumulated volume in the coolant hose (increased dripping possible)

Observe hygiene to prevent infection: The liquid in the coolant hose may not draw back to far (avoid contamination).

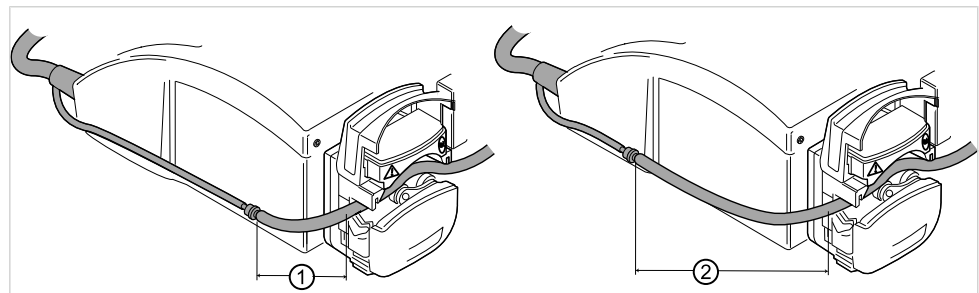

7.2 Exchange coolant hose




Note

The parts delivered for exchange are not sterile. Rise and sterilise them before the first treatment.

The coolant hose is run to the handpiece/contra-angle handpiece and consists of 2 pieces:


Coolant hose ① (approx. 2.3 m long and obtainable as a section under **Mat. no. 0.065.5279**)

Coolant hose ④ (obtainable approx. 150 mm long under **Mat. no. 0.593.0252** and as a section under **Mat. no. 0.065.5188**)

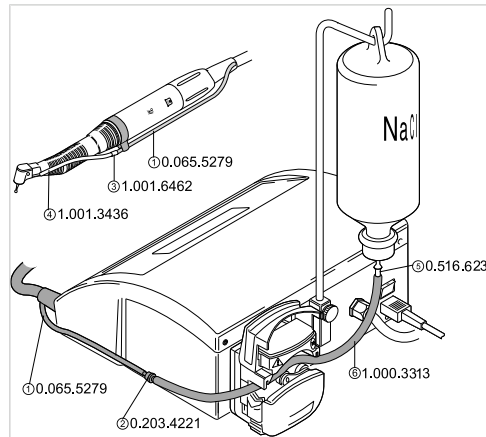

7.2.1 Exchange the coolant hose from the handpiece or angle piece.

- ▶ Pull off the coolant hose ④ **Mat. no. 0.593.0252** (150 mm), or **Mat. no. 1.001.3436** (70 mm), or **Mat. no. 0.065.5188** (cut section) off the coupling ③ **Mat. no. 1.001.6462** and the handpiece or angle piece and exchange.

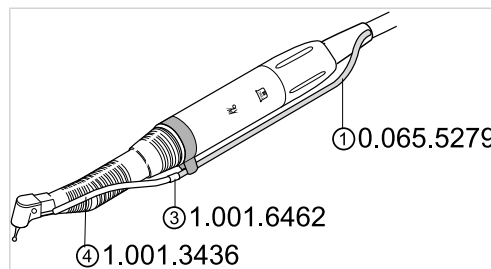

7.2.2 Exchange the coolant hose in the motor hose




Note


Run the new hose straight in the direction of tension so that it can be pulled in as easily as possible. Talcum powder can be used as a lubricant.

- ▶ Remove the coolant hose ① from the plug-in nipple ② and coupling ③.

- ▶ Connect the new and existing coolant hose ① to the coupling ③.

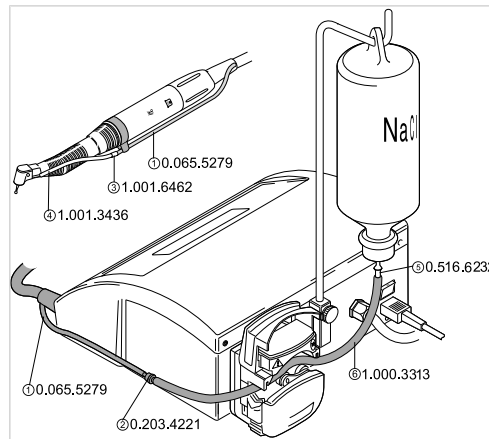

- ▶ Carefully remove the old hose from the motor end and simultaneously shove on the new one until the coupling ③ piece is visible.
- ▶ Pull the coupling ③ out of the hose opening e.g. using a tweezers so that the coupling piece is not clamped there.

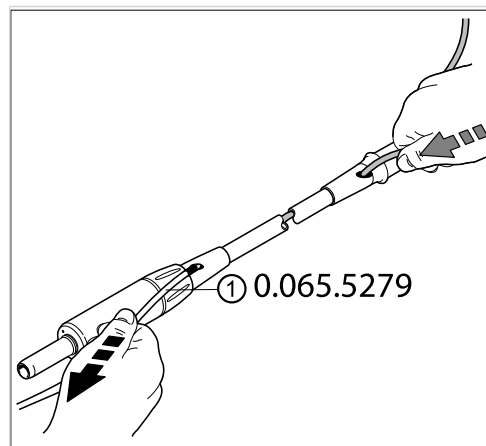

7.2.3 Disposable coolant hose


Sterile-packaged coolant hoses can be ordered under **Mat. no. 1.001.9902** (package contents: 10 units).

This should be used for applications with an increased infection risk or when the motor hose cannot be prepared according to the existing instructions.

The disposable coolant hose consists of an insertion tip, pump hose, clip and hose holder.

The length of the hose can be cut at the exit end with a sterile knife or sterile scissors. Observe the instructions regarding use on the package.

- ▶ Connect the coolant hose directly to the outside of the motor hose without loops or kinks, and affix with the accompanying holders (evenly distributed).

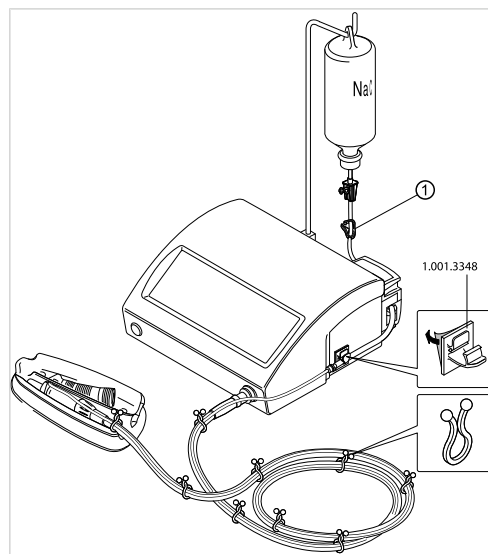

- ▶ Before using the pump, completely open the inserted hose clip ①.

The coolant flow is also affected by the outlet in the handpiece or angle piece. A focused water jet should arise.

- ▶ Set the pump flow on the device.



Note


To avoid dripping, keep the pump hose length after the pump as short as possible.

7.3 Device settings table

Enter the set parameters in the table before each service check. Always check the parameters using the table after the service check.

Date:		Name:			
Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
Free					
Date:		Name:			
Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
Free					
Date:		Name:			
Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
Free					
Date:		Name:			
Step	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
	Instrument	Speed (rpm)	Torque (Ncm)	Pump	Motor r/l
Free					

8 Accessories

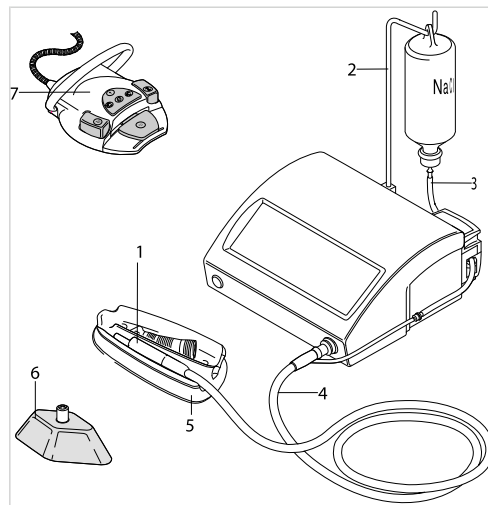

The INTRAsurg 300/300 plus contains the follow permissible accessories:

- INTRA S 550 surgical motor **Mat. no. 1.000.8072** ① (autoclavable collectorless motor with a speed range of 300 to 40000 rpm)
- Alternative: INTRA LUX SL 550 surgical motor **Mat. no. 1.001.3421** ① (autoclavable collectorless motor with a speed range of 300 to 40000 rpm)

- Bottle holder **Mat. no. 0.761.1872** ②
Do not use any container with a volume in excess of 1 litre (danger of tipping over).
- Pump hose **Mat. no. 1.001.3347** ③
- Disposable coolant hose sterile **Mat. no. 1.001.9902**
- Motor hose 2 m (standard equipment) **Mat. no. 1.001.2651** ④
Motor hose 3 m (please inquire) **Mat. no. 1.004.6825** ④
The connection hose with coolant guidance can be autoclaved.
- Motor-Steri support **Mat. no. 0.726.2922** ⑤
- IPX 8 foot control element with 3.5 m connection cable (please inquire) **Mat. no. 1.004.8276** ⑦
- Steri-Set INTRAsurg 300 (please inquire) **Mat. no. 1.001.4953** includes the following parts: 1x motor hose 2 m **Mat. no. 1.001.2651**; 1x surgical motor INTRA S 550; 1x motor-Steri support **Mat. no. 0.726.2922**; 1x hose 150mm **Mat. no. 0.593.0252**; 1x pump hose **Mat. no. 1.001.3347**.

For the INTRAsurg 300 Plus:

- INTRAsurg Calibration **Mat. no. 1.002.3569** ⑥ (disinfectable).
- Steri-Set INTRAsurg 300plus (only upon request) **Mat. no. 1.001.4968** consisting of the following parts: 1 motor hose 2 m **Mat. no. 1.001.2651**; 1 surgical motor INTRA LUX SL 550; 1 motor sterile holder **Mat. no. 0.726.2922**; 1 hose 150 mm **Mat. no. 0.593.0252**; 1 pump hose **Mat. no. 1.001.3347**.

KaVo surgical, handpiece and angle piece line corresponding to the order.
For news on surgery and implantology, see www.kavo.com.

9 Details on electromagnetic compatibility

9.1 Guidelines and manufacturer's declaration - electromagnetic transmission

The IS 300 is for use in an environment like the one cited below. The IS 300 customer or user should ensure that use takes place in such an environment.

Measurements of noise transmissions	Conformance	Electromagnetic environment - hints
HF-transmissions according to CISPR 11	Group 1	The IS 300 uses HF energy only for its internal operation. Its HF transmission is therefore very low, and it is improbable that neighbouring electronic devices will be disturbed.
HF-transmissions according to CISPR 11	Class B	The IS 300 is for use in all facilities including residential ones, and facilities that are directly connected to a public power supply that also supplies residential buildings.
Transmissions of harmonics according to IEC 61000-3-2	Class A	The IS 300 is for use in all facilities including residential ones, and facilities that are directly connected to a public power supply that also supplies residential buildings.
Transmissions of voltage fluctuations or flicker according to IEC 61000-3-3	fulfilled	The IS 300 is for use in all facilities including residential ones, and facilities that are directly connected to a public power supply that also supplies residential buildings.



Note

The device or system may not be used or stacked directly next to other devices. If it has to be used close to or stacked next to other devices, the device or system must be monitored to ensure that it is used properly in the existing arrangement.



Note

The immunity test levels required in IEC 60601 are met.

9.2 Guidelines and manufacturer's declaration - electromagnetic resistance to jamming


The IS 300 is for use in an environment like the one cited below. The IS 300 customer or user should ensure that use takes place in such an environment.

Immunity tests	IEC 60601 test level	Conformance level	Electromagnetic environment - guidelines
Electrostatic discharge (ESD) according to IEC 61000-4-2	± 6 kV contact discharge ± 8 kV atmospheric discharge	± 6 kV contact discharge ± 8 kV atmospheric discharge	Floors should be made of wood or concrete or have ceramic tiles. When the floor is made of synthetic material, the relative humidity must be at least 30%.
Fast transient electrical disturbances/ Bursts according to IEC 61000-4-4	± 2 kV for power lines ± 1 kV for input and output lines	± 2 kV for power lines ± 1 kV for input and output lines	The quality of the supply voltage should correspond to that of a typical business or hospital environment.
Surges according to IEC 61000-4-5	± 1 kV Push-pull voltage (symmetrical) ± 2 kV common mode voltage (unsymmetrical)	± 1 kV Push-pull voltage (symmetrical) ± 2 kV common mode voltage (unsymmetrical)	The quality of the supply voltage should correspond to that of a typical business or hospital environment.
Voltage interruptions, short-term interruptions and fluctuations of the supply voltage according to IEC 61000-4-11	< 5% U_T for 1/2 period (>95% interruption) 40 % U_T for 5 periods (60% interruption) 70 % U_T for 25 periods (30 % interruption) < 5% U_T for 5 s (>95% interruption)	< 5% U_T for 1/2 period (>95% interruption) 40 % U_T for 5 periods (60% interruption) 70 % U_T for 25 periods (30 % interruption) < 5% U_T for 5 s (>95% interruption)	The quality of the supply voltage should correspond to that of a typical business or hospital environment. When the user of the IS 300 needs continued operation even when the power supply is interrupted, it is recommended to supply the IS 300 from an uninterrupted power supply or a battery.
Magnetic field with a supply frequency (50/60 Hz) according to IEC 61000-4-8	3 A/m	3 A/m	Magnetic fields at the mains frequency should correspond to typical values in a business and hospital environment.

NOTE: U_T is the alternating mains voltage before the test level is used.

9.3 Guidelines and manufacturer's declaration - electromagnetic resistance to jamming

The IS 300 is for use in an environment like the one cited below. The IS 300 customer or user should ensure that use takes place in such an environment.

Immunity tests	IEC 60601 test level	Conformance level	Electromagnetic environment - guidelines
Conducted HF disturbances according to IEC 61000-4-6 Radiated HF disturbances according to IEC 61000-4-3	3 V _{eff} 150 kHz to 80 MHz 10 V/m 80 MHz to 2.5 GHz	3 V _{eff} 3 V/m	Portable and mobile radio devices should not be used closer to the IS 300 including the wires, than the recommended safe distance calculated using the equation for the transmission frequency. Recommended safe distance: $d = [3.5/3]\sqrt{P} = 1.17\sqrt{P}$ $d = [3.5/3]\sqrt{P} = 1.17\sqrt{P}$ for 80 MHz to 800 MHz $d = [7.0/3]\sqrt{P} = 2.33\sqrt{P}$ for 800 MHz to 2.5 GHz with P as the maximum rated power of the transmitter in Watts (W) according to the transmitter manufacturer, and d as the recommended safe distance in meters (m). The field strength of stationary radio transmitters should be less than the conformance level at all frequencies in an on-site check ^{a, b} Disturbances are possible close to devices that have the following symbol. 

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not be applicable in every case. The spread of electromagnetic waves is absorbed and reflected by buildings, objects and people.

^aThe field strength of stationary transmitters such as base stations of mobile telephones and land radio devices, amateur radio stations, AM and FM, radio and television broadcasters cannot be theoretically predetermined. To determine the electromagnetic environment of stationary transmitters, a study of the location should be considered. When the measured field strength at the location on which the IS 300 is used, exceeds the conformity level, the IS 300 should be watched to ensure that it is functioning as per the correct usage. Should unusual performance features be observed, additional measures may be required, such as e.g. a different alignment or another location for the IS 300.

^bWithin the frequency range of 150 kHz to 80 MHz, the field strength should be less than 3V_{eff} V/m.

9.4 Recommended safe distance between portable and mobile HF telecommunications equipment and the INTRAsurg 300

The IS 300 is for use in an environment like the one cited below. The customer or user of the IS 300 can help to avoid electromagnetic faults by keeping to the minimum distance apart between portable and mobile HF-telecommunication devices (transmitters) and the IS 300 - dependent on the output lines for the communication device - as given below.

Rated power of the transmitter in W	150 kHz to 80 MHz $d=1.17\sqrt{P}$	80 MHz to 800 MHz $d=1.17\sqrt{P}$	800 MHz to 2.5 GHz $d=2.33\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,37	0,37	0,74
1	1,17	1,17	2,33
10	3,70	3,70	7,37
100	11,70	11,70	23,30

For transmitters whose maximum rated power is not in the above table, the recommended safe distance d in meters (m) can be calculated using the equation for the respective gap, where P is the maximum rated power of the transmitter in Watts (W) according to the manufacturer's information.

NOTE 2: These guidelines may not be applicable in every case. The spread of electromagnetic waves is absorbed and reflected by buildings, objects and people.

Comment 1: To calculate the recommended safe distance from transmitters with a frequency range of 80 MHz to 2.5 GHz, an additional factor of 10/3 was used to reduce the probability that a mobile unintentionally brought into the patient area would cause malfunction.



KaVo. Dental Excellence.